

IUI MEMO

UTERINE INSEMINATOR FOR WASHED SPERMATOZOA WITH MEMORY FUNCTION

The IUI Memo is an intrauterine insemination cannula with MEMO function, which proactively enables the distal tip to bend according to cervical flexion.

FEATURES


- positioner to prevent damage to fundus
- two lateral ports at distal end for optimal sperm distribution
- arrows on handle to indicate lateral ports at distal end
- markings on outer tube
- memory function
- handle with luer lock for a more secure fitting to syringes



ORIGINAL PICTURE

GYNO INFO IUI MEMO

Length	206 mm (± 2 mm)
Tube OD	1.85 mm (± 0.10 mm)
Sterile	sterilized by gamma irradiation
Packaging	125/shipping box
Quality control	MEA and LAL
Order number	#4220 Memo 

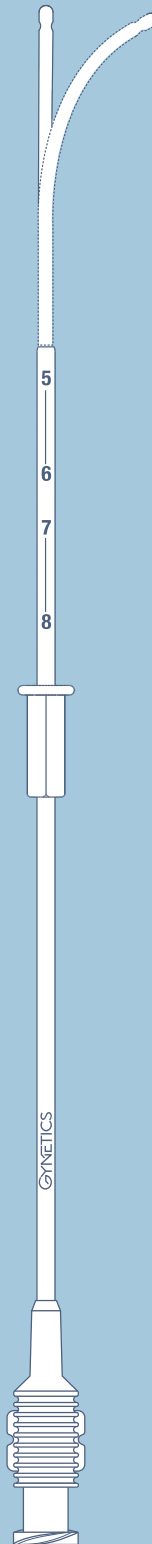

GMP Belgium IUI Memo 09/2014



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine