

THE TULIP MEMO

EMBRYO TRANSFER CATHETER WITH TULIP SHAPED BULB

The Tulip is a set for difficult cervical passage, consisting of two parts: a guiding catheter (A) and a loading catheter (C).


FEATURES


- the Tulip tip enables easy atraumatic entry through the cervical canal
- mandrel available for extra difficult transfers
- positioner for depth adjustment
- a combination of a soft, flexible intrauterine catheter and a solid cervix catheter with memory
- markings on outer tube for easy positioning
- memo function of the guiding catheter
- 5 free mandrels to preshape the guiding catheter in every box of 25 pieces
- handle with luer lock fitting for a more secure connection to syringes
- reinforced metal support of loading catheter for optimal control



ORIGINAL PICTURE

GYNO INFO THE TULIP MEMO

Length	A) 200 mm (\pm 2 mm) C) 266 mm (\pm 2 mm)
Tube OD	A) 2.20 mm (\pm 0.05 mm) C) 1.00 mm (\pm 0.05 mm)
Total volume	C) 0.21 ml
Total volume minus hub	C) 0.069 ml
Tulip bulb	3 mm
Sterile	sterilized by radiation
Packaging	75/shipping box
Quality control	MEA and LAL
Order number	#4000 Memo 

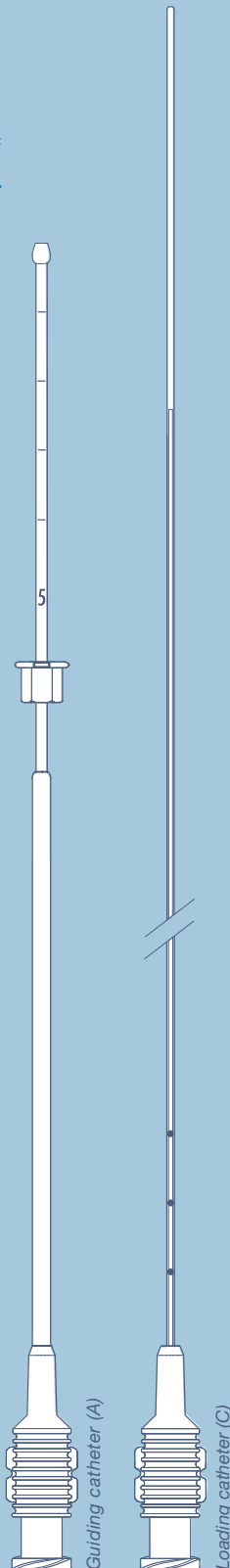

GMP Belgium The Tulip Memo 04/2011