


GYNETICS®

Practise better Medicine

RENEWED EMBRYO TRANSFER CATHETER PORTFOLIO

The renewed Embryo Transfer catheters offer you a variety of possibilities to enhance your practice

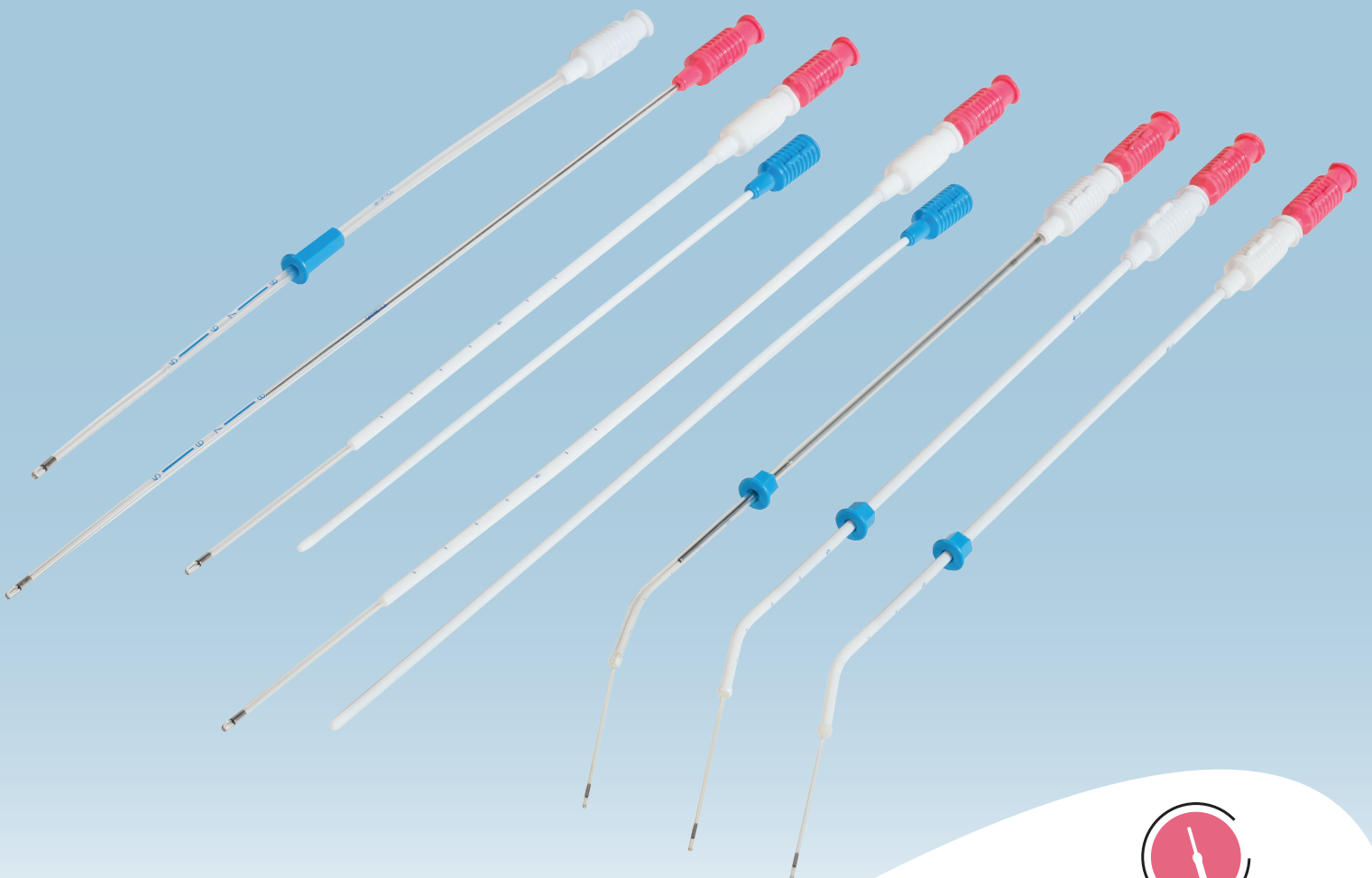

Embryo Transfer


The Tulip

EMBRYO TRANSFER CATHETER SETS WITH BULB SHAPED TIP


The renewed Tulip Embryo Transfer catheter range offers a variety of possibilities to support you enhancing your practice.

- The **bulb shaped** guiding catheters assist atraumatic entry into the cervical canal and are available in different shapes.
- The **Memo function** allows shaping the guiding catheter in any desired silhouette to suit various patients' anatomies.
- The Tulip **loading catheters** are available with a Classic tip or a Tip-in-Sight (TIS), which supports precise positioning under ultrasound guidance.

Guiding Catheter **bulb shape** options

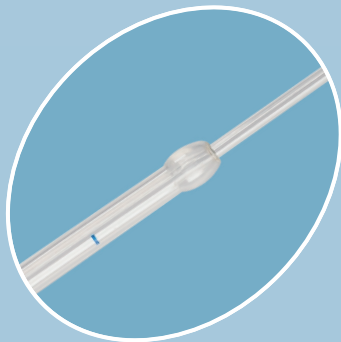

The Tulip S
Oval shaped bulb

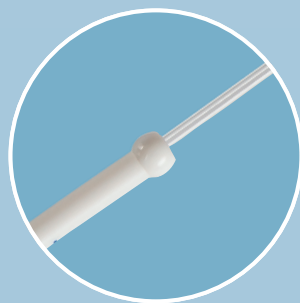


The Tulip L
Rigid round shaped bulb



The Tulip LS
SOFT round shaped bulb

Loading Catheter **tip** options

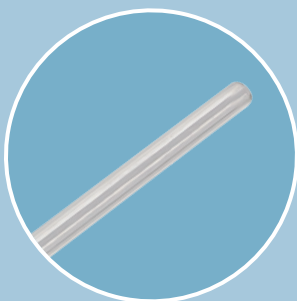

Classic




Tip-in-Sight (TIS)
Echogenic

CHOOSE YOUR TULIP SET

PRODUCT	ORDER NUMBER	FEATURES				
		Echogenic tip	Preshaped with curve indicator	Memory function	Mandrel available 5 per 25 sets	
 <p>Oval shaped bulb</p>	THE TULIP S	#4000 Tulip S		✓		
	THE TULIP S Memo	#4000 S Memo			✓	✓
	THE TULIP S Tip-in-Sight	#4000 Tulip S TIS	✓	✓		
	THE TULIP S Memo Tip-in-Sight	#4000 S Memo TIS	✓		✓	✓
 <p>Rigid round shaped bulb</p>	THE TULIP L	#4000 Tulip L		✓		
	THE TULIP L Memo	#4000 L Memo			✓	✓
	THE TULIP L Tip-in-Sight	#4000 Tulip L TIS	✓	✓		
	THE TULIP L Memo Tip-in-Sight	#4000 L Memo TIS	✓		✓	✓
 <p>SOFT round shaped bulb</p>	THE TULIP LS	#4000 Tulip LS		✓		
	THE TULIP LS Memo	#4000 LS Memo			✓	✓
	THE TULIP LS Tip-in-Sight	#4000 Tulip LS TIS	✓	✓		
	THE TULIP LS Memo Tip-in-Sight	#4000 LS Memo TIS	✓		✓	✓

GYNO INFO

Length	guiding catheter: 200 mm loading catheter: 266 mm
Tube OD	guiding catheter: 2.20 mm loading catheter: 1.00 mm
Total volume Total volume minus hub	loading catheter: 0.21 ml loading catheter: 0.069 ml
Bulb shape	oval shaped: 3 mm rigid/soft round shaped: 2.3 mm
Sterile	sterilized by gamma irradiation
Packaging	75/shipping box
Quality control	MEA and LAL

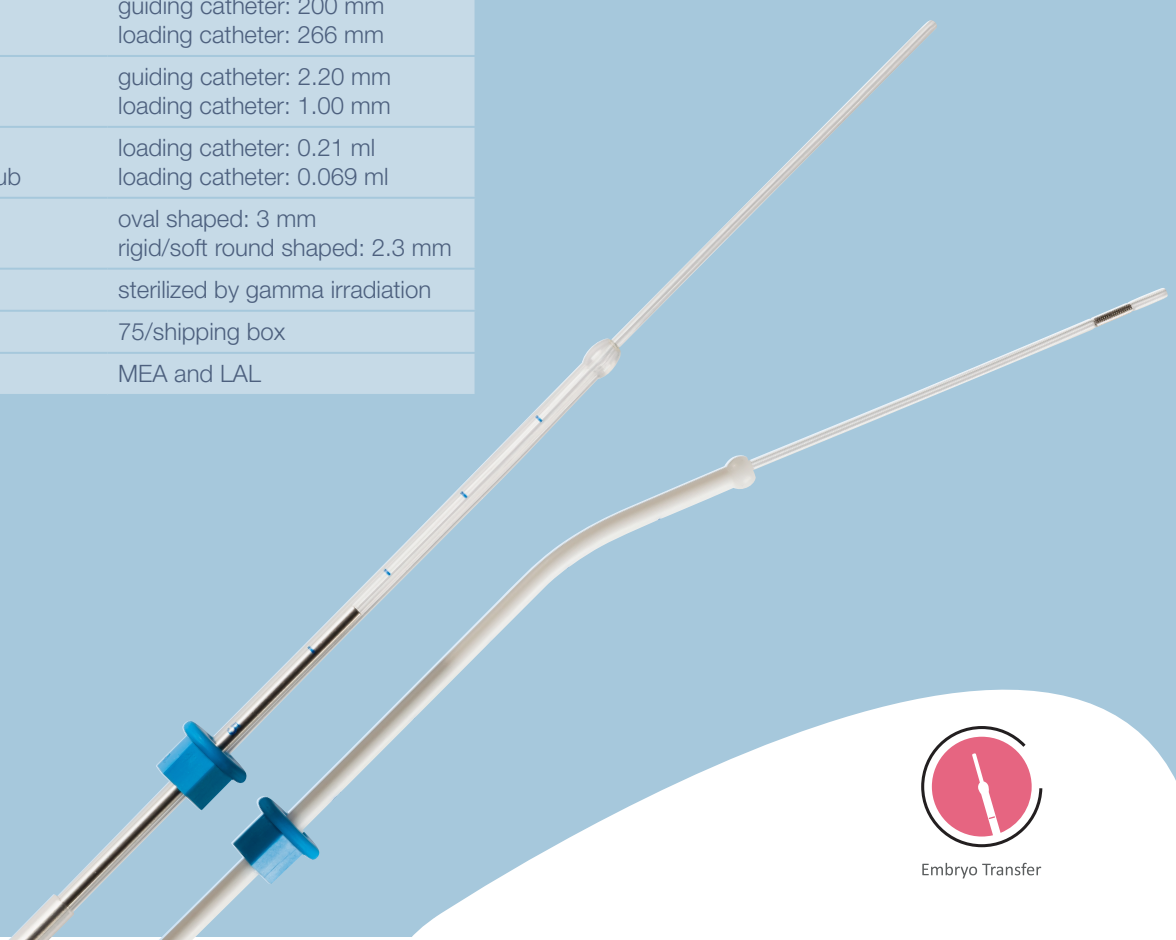


Embryo Transfer

Emtrac Delphin / Semtrac / The Tulip

EMBRYO TRANSFER CATHETERS WITH TIP-IN-SIGHT (TIS) TECHNOLOGY

The TIS technology supports precise positioning of the loading catheter tip under ultrasound guidance. The new range of Embryo Transfer catheters with the TIS technology is an extension of the existing portfolio with Classic loading tips.

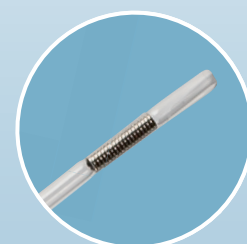

GYNO INFO

The Tulip TIS EMBRYO TRANSFER CATHETER SETS WITH BULB SHAPED TIP								
PRODUCT	ORDER NUMBER	GUIDING CATHETER		LOADING CATHETER		VOLUME		
		LENGTH	TUBE OD	LENGTH	TUBE OD	TOTAL	TOTAL MINUS HUB	BULB SHAPE
THE TULIP S Tip-in-Sight	#4000 Tulip S TIS	200 mm	2.20 mm	266 mm	1.00 mm	0.21 ml	0.069 ml	3 mm
THE TULIP S Memo Tip-in-Sight	#4000 S Memo TIS	200 mm	2.20 mm	266 mm	1.00 mm	0.21 ml	0.069 ml	3 mm
THE TULIP L Tip-in-Sight	#4000 Tulip L TIS	200 mm	2.20 mm	266 mm	1.00 mm	0.21 ml	0.069 ml	2.3 mm
THE TULIP L Memo Tip-in-Sight	#4000 L Memo TIS	200 mm	2.20 mm	266 mm	1.00 mm	0.21 ml	0.069 ml	2.3 mm
THE TULIP LS Tip-in-Sight	#4000 Tulip LS TIS	200 mm	2.20 mm	266 mm	1.00 mm	0.21 ml	0.069 ml	2.3 mm
THE TULIP LS Memo Tip-in-Sight	#4000 LS Memo TIS	200 mm	2.20 mm	266 mm	1.00 mm	0.21 ml	0.069 ml	2.3 mm

Semtrac TIS SOFT EMBRYO TRANSFER CATHETER								
PRODUCT	ORDER NUMBER	GUIDING CATHETER		LOADING CATHETER		VOLUME		
		LENGTH	TUBE OD	LENGTH	TUBE OD	TOTAL	TOTAL MINUS HUB	
Semtrac C Tip-in-Sight	#2000 C TIS	-	-	252 mm	1.70 mm	0.20 ml	0.065 ml	
Semtrac Set 4.2 Tip-in-Sight	#2000 Set 4.2 TIS	158 mm	2.55 mm	232 mm	1.70 mm	0.23 ml	0.081 ml	
Semtrac Set 5 Tip-in-Sight	#2000 Set 5 TIS	178 mm	2.55 mm	252 mm	1.70 mm	0.22 ml	0.089 ml	
Semtrac Set W-180 Tip-in-Sight	#2000 Set W-180 TIS	141 mm	2.55 mm	214 mm	1.70 mm	0.23 ml	0.074 ml	
Semtrac Set W-240 Tip-in-Sight	#2000 Set W-240 TIS	203 mm	2.55 mm	277 mm	1.70 mm	0.26 ml	0.098 ml	

Emtrac Delphin TIS EMBRYO TRANSFER CATHETER WITH POSITIONER								
PRODUCT	ORDER NUMBER	GUIDING CATHETER		LOADING CATHETER		VOLUME		
		LENGTH	TUBE OD	LENGTH	TUBE OD	TOTAL	TOTAL MINUS HUB	
Emtrac Delphin Tip-in-Sight	#4219 Delphin TIS	-	-	210 mm	2.00 mm	0.23 ml	0.072 ml	

BEFORE USING REFER TO INSTRUCTIONS FOR USE (IFU) FOR INDICATIONS, CONTRA-INDICATIONS AS WELL AS WARNINGS AND PRECAUTIONS.

90240118_001