

GYNO INFO

INTRAUTERINE INSEMINATION CANNULA RANGE

PRODUCT	ORDER NUMBER
---------	--------------

IUI	#4220
IUI Memo	#4220 Memo
IUI Plus	#4220 Plus
Smooze IUI	#4225
Smooze IUI Medium	#4225 M
Smooze IUI Long	#4225 L

GMP Belgium IUI Range 06/2019



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine

IUI

UTERINE INSEMINATOR FOR WASHED SPERMATOZOA

The IUI is an intrauterine insemination cannula with two lateral distal ports for optimal sperm distribution.

FEATURES

- positioner to prevent damage to fundus
- two lateral ports at distal end for optimal sperm distribution
- arrows on handle to indicate lateral ports at distal end
- markings on the outer tube
- handle with luer lock for a more secure fitting to syringes



ORIGINAL PICTURE

GYNO INFO IUI

Length	206 mm
Tube OD	1.85 mm
Sterile	sterilized by gamma irradiation
Packaging	125/shipping box
Quality control	MEA and LAL
Order number	#4220



GMP Belgium IUI 09/2014



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine

IUI MEMO

UTERINE INSEMINATOR FOR WASHED SPERMATOZOA WITH MEMORY FUNCTION

The IUI Memo is an intrauterine insemination cannula with MEMO function, which proactively enables the distal tip to bend according to cervical flexion.

FEATURES


- positioner to prevent damage to fundus
- two lateral ports at distal end for optimal sperm distribution
- arrows on handle to indicate lateral ports at distal end
- markings on outer tube
- memory function
- handle with luer lock for a more secure fitting to syringes



ORIGINAL PICTURE

GYNO INFO IUI MEMO

Length	206 mm
Tube OD	1.85 mm
Sterile	sterilized by gamma irradiation
Packaging	125/shipping box
Quality control	MEA and LAL
Order number	#4220 Memo

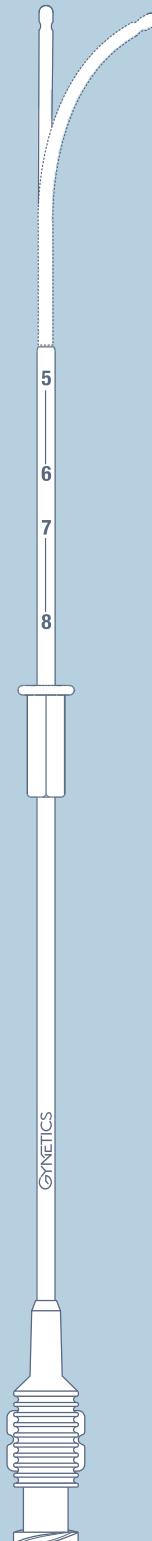

GMP Belgium IUI Memo 09/2014



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine

IUI PLUS

UTERINE INSEMINATOR FOR WASHED SPERMATOZOA WITH MANDREL

The IUI Plus is a set for difficult cervical passage, consisting of two parts:

- A) Loading cannula
- B) Mandrel


FEATURES

- positioner for depth adjustment
- a combination of a soft, flexible intrauterine cannula and a solid cervix cannula
- markings on outer tube for easy positioning
- mandrel to prevent mucus to enter the loading cannula
- handle with luer lock for a more secure fitting to syringes

ORIGINAL PICTURE

GYNO INFO IUI PLUS

Length	A) 216 mm B) 240 mm
Tube OD	A) 1.60 mm B) 1.00 mm
Total volume	A) 0.42 ml
Total volume minus hub	A) 0.258 ml
Sterile	sterilized by gamma irradiation
Packaging	75/shipping box
Quality control	MEA and LAL
Order number	#4220 Plus

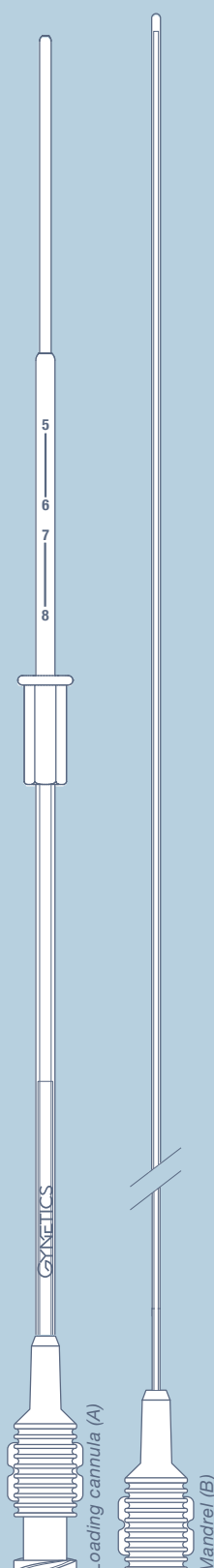

GMP Belgium IUI Plus 06/2014



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS®
Practise better Medicine

SMOOZE IUI

FLEXIBLE INTRAUTERINE INSEMINATION CANNULA

The Smooze IUI is a flexible intrauterine insemination cannula with a smooth tip for less traumatic inseminations.

FEATURES

- smooth tip with one distal side hole
- very flexible and therefore less traumatic
- handle with luer lock for a more secure fitting to syringes
- available in 3 lengths:


- A) Short 7 cm (not including handle)
B) Medium 14 cm (not including handle)
C) Long 20 cm (not including handle)



ORIGINAL PICTURE

GYNO INFO SMOOZE IUI

Length :	A) 100 mm B) 169 mm C) 230 mm
Tube OD :	2 mm
Volume :	A) min. 0.030 ml B) min. 0.057 ml C) min. 0.080 ml
Sterile :	ETO sterilized
Packaging :	A) 125/shipping box B) 125/shipping box C) 75/shipping box
Quality control :	MEA and LAL
Order number :	A) #4225 B) #4225 M C) #4225 L

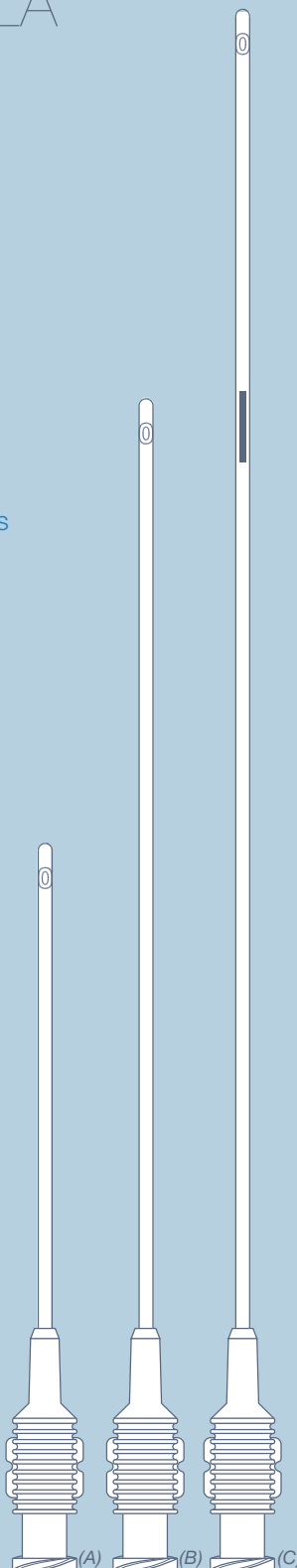

GMP Belgium Smooze IUI 09/2014



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine

GYNO INFO

EMBRYO TRANSFER CATHETER RANGE

EMBRYO TRANSFER CATHETERS FOR EASY TRANSFER

PRODUCT	ORDER NUMBER
Emtrac A	#4219 A
Emtrac Delphin	#4219 Delphin
Emtrac Delphin Tip-in-Sight	#4219 Delphin TIS
Semtrac C	#2000 C
Semtrac C Tip-in-Sight	#2000 C TIS

EMBRYO TRANSFER CATHETER SETS FOR DIFFICULT TRANSFER

PRODUCT	ORDER NUMBER
Emtrac Set	#4219 Set
Semtrac Set 4.2	#2000 Set 4.2
Semtrac Set 4.2 Tip-in-Sight	#2000 Set 4.2 TIS
Semtrac Set 5	#2000 Set 5
Semtrac Set 5 Trial	#2000 Set 5 T
Semtrac Set 4.2 Tip-in-Sight	#2000 Set 4.2 TIS
Semtrac Set W-180	#2000 W-180
Semtrac Set W-180 Trial	#2000 W-180 T
Semtrac Set W-180 Tip-in-Sight	#2000 W-180 TIS
Semtrac Set W-240	#2000 W-240
Semtrac Set W-240 Trial	#2000 W-240 T
Semtrac Set W-180 Tip-in-Sight	#2000 W-180 TIS

GMP Belgium ET Range 06/2019



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine

GYNO INFO

EMBRYO TRANSFER CATHETER RANGE

EMBRYO TRANSFER CATHETER SETS FOR DIFFICULT TRANSFER

PRODUCT	ORDER NUMBER
The Tulip	#4000 The Tulip
The Tulip Memo	#4000 Memo
The Tulip S	#4000 Tulip S
The Tulip S Memo	#4000 Tulip S Memo
The Tulip S Tip-in-Sight	#4000 Tulip S TIS
The Tulip S Memo Tip-in-Sight	#4000 Tulip S Memo TIS
The Tulip L	#4000 Tulip L
The Tulip L Memo	#4000 Tulip L Memo
The Tulip L Tip-in-Sight	#4000 Tulip L TIS
The Tulip L Memo Tip-in-Sight	#4000 Tulip L Memo TIS
The Tulip LS	#4000 Tulip LS
The Tulip LS Memo	#4000 Tulip LS Memo
The Tulip LS Tip-in-Sight	#4000 Tulip LS TIS
The Tulip LS Memo Tip-in-Sight	#4000 Tulip LS Memo TIS

GMP Belgium ET Range 06/2019



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine

EMTRAC A

EMBRYO TRANSFER CATHETER WITH POSITIONER

The Emtrac A is a combination of a soft, flexible intrauterine catheter and a solid cervix catheter.


FEATURES

- positioner for depth adjustment
- a combination of a soft, flexible intrauterine catheter and a solid cervix catheter
- markings on outer tube for easy positioning
- handle with luer lock for a more secure fitting to syringes

ORIGINAL PICTURE

GYNO INFO EMTRAC A

Length	216 mm
Tube OD	1.60 mm
Total volume	0.42 ml
Total volume minus hub	0.258 ml
Sterile	sterilized by radiation
Packaging	125/shipping box
Quality control	MEA and LAL
Order number	#4219 A

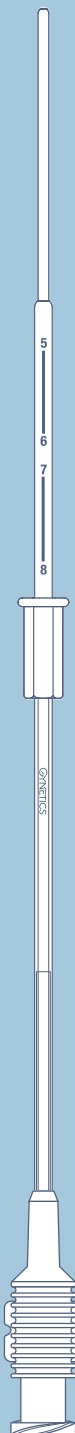

GMP Belgium Emtrac A 04/2011



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine

EMTRAC DELPHIN

EMBRYO TRANSFER CATHETER WITH POSITIONER AND SOFT UTERINE CATHETER

The combination of a soft flexible intrauterine catheter and a solid cervix catheter creates an universal catheter for easy insertion in the most difficult cervical tract.

FEATURES

- positioner to prevent damage to fundus
- one piece catheter
- markings to indicate depth
- extra soft uterine catheter to ensure safe insertion through cervix
- handle with luer lock for a more secure fitting to syringes
- **available with Tip-in-Sight (TIS) technology for precise positioning of the catheter tip under ultrasound guidance**

TIP OPTIONS



GYNO INFO EMTRAC DELPHIN

Length	210 mm
Tube OD	2.00 mm
Total volume	0.23 ml
Total volume minus hub	0.072 ml
Sterile	ETO sterilized
Packaging	75/shipping box
Quality control	MEA and LAL
Order numbers	
- Emtrac Delphin	#4219 Delphin
- Emtrac Delphin Tip-in-Sight	#4219 Delphin TIS

BEFORE USING REFER TO INSTRUCTIONS FOR USE (IFU) FOR INDICATIONS, CONTRA-INDICATIONS AS WELL AS WARNINGS AND PRECAUTIONS.



SEMTRAC C

SOFT EMBRYO TRANSFER CATHETER

The Semtrac C is a soft embryo transfer catheter.

FEATURES

- soft bullet shaped tip to minimize trauma
- one piece economical catheter
- markings to indicate depth
- metal tube inside for better control
- handle with luer lock for a more secure fitting to syringes
- **available with Tip-in-Sight (TIS) technology for precise positioning of the catheter tip under ultrasound guidance**

TIP OPTIONS



CLASSIC

TIP-IN-SIGHT

GYNO INFO SEMTRAC C

Length	252 mm
Tube OD	1.70 mm
Total volume	0.20 ml
Total volume minus hub	0.065 ml
Sterile	ETO sterilized
Packaging	75/shipping box
Quality control	MEA and LAL
Order numbers	
- Semtrac C	#2000 C
- Semtrac C Tip-in-Sight	#2000 C TIS

BEFORE USING REFER TO INSTRUCTIONS FOR USE (IFU) FOR INDICATIONS, CONTRA-INDICATIONS AS WELL AS WARNINGS AND PRECAUTIONS.



EMTRAC SET

EMBRYO TRANSFER CATHETER WITH POSITIONER

The Emtrac Set is a set for difficult cervical passage, consisting of three parts:

- A) Guiding catheter
- B) Mandrel
- C) Loading catheter


FEATURES

- positioner for depth adjustment
- a combination of a soft, flexible intrauterine catheter and a solid cervix catheter
- markings on outer tube for easy positioning
- soft loading catheter to minimize trauma
- handle with luer lock for a more secure fitting to syringes

ORIGINAL PICTURE

GYNO INFO EMTRAC SET

Length	A) 216 mm B) 240 mm C) 266 mm
Tube OD	A) 1.60 mm B) 1.00 mm C) 1.00 mm
Total volume	C) 0.20 ml
Total volume minus hub	C) 0.069 ml
Sterile	sterilized by radiation
Packaging	75/shipping box
Quality control	MEA and LAL
Order number	#4219 Set

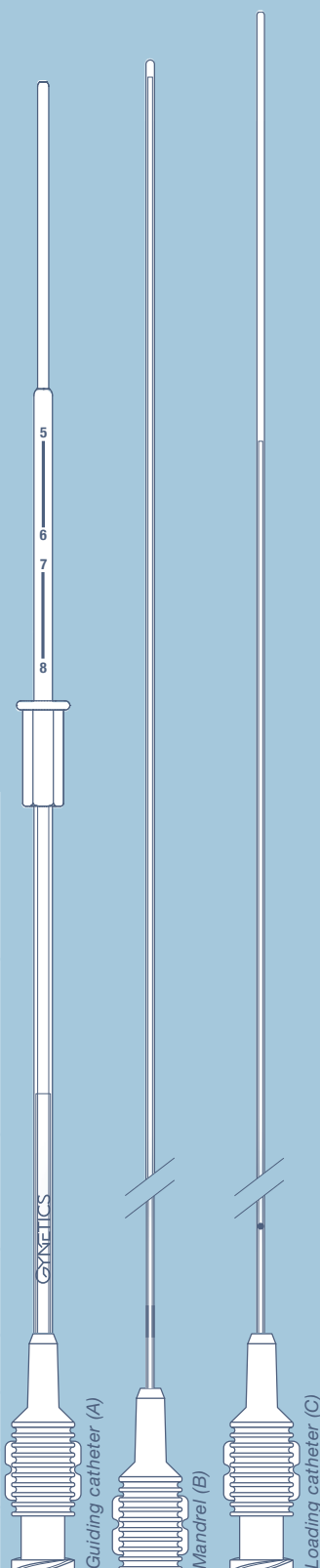

GMP Belgium Emtrac Set 04/2011



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine

SEMTRAC SET 4.2

SOFT EMBRYO TRANSFER CATHETER

The Semtrac Set is a soft embryo transfer catheter set for difficult embryo transfers.

FEATURES

- soft bullet shaped tip to minimize trauma
- three piece catheter for difficult transfers
- 5 **free** mandrels (B), to preshape the guiding catheter, in every box of 25 pieces
- reinforced shaft of loading catheter for optimal control
- handle with luer lock for a more secure fitting to syringes
- **available with Tip-in-Sight (TIS) technology for precise positioning of the catheter tip under ultrasound guidance**


LOADING CATHETER TIP OPTIONS



GYNO INFO SEMTRAC SET 4.2

Guiding catheter (A) length	158 mm
tube OD	2.55 mm
Loading catheter (C) length	232 mm
tube OD	1.70 mm
total volume	0.23 ml
total volume minus hub	0.081 ml
Sterile	ETO sterilized
Packaging	75/shipping box
Quality control	MEA and LAL
Order numbers	
- Semtrac Set 4.2	#2000 Set 4.2
- Semtrac Set 4.2 Tip-in-Sight	#2000 Set 4.2 TIS

BEFORE USING REFER TO INSTRUCTIONS FOR USE (IFU) FOR INDICATIONS, CONTRA-INDICATIONS AS WELL AS WARNINGS AND PRECAUTIONS.

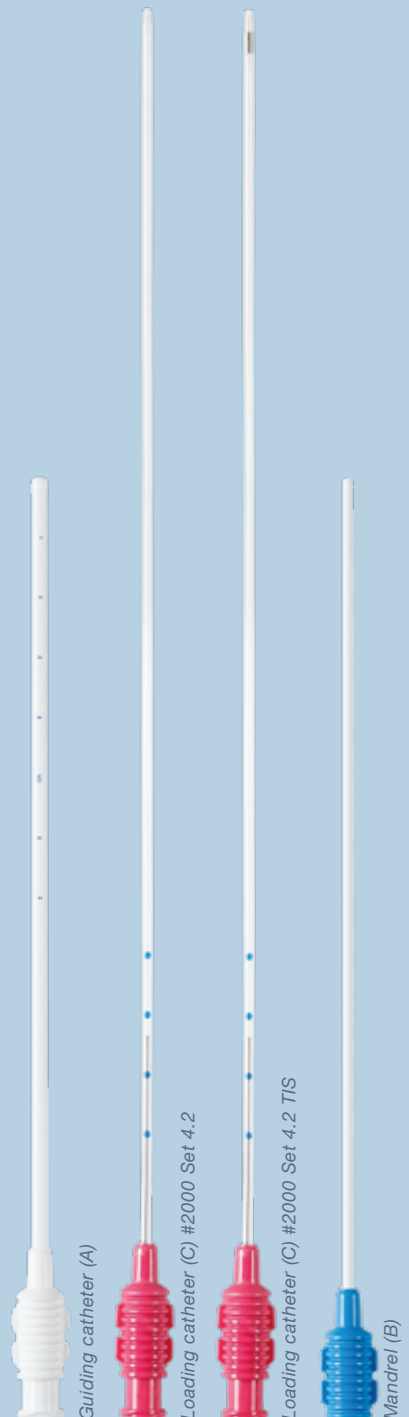

SEMTRAC SET 5

SOFT EMBRYO TRANSFER CATHETER

The Semtrac Set is a soft embryo transfer catheter set for difficult embryo transfers.

FEATURES

- soft bullet shaped tip to minimize trauma
- three piece catheter for difficult transfers
- 5 **free** mandrels (B), to preshape the guiding catheter, in every box of 25 pieces
- reinforced shaft of loading catheter for optimal control
- handle with luer lock for a more secure fitting to syringes
- trial catheter set available for cervix examination prior to embryo transfer consisting of a guiding catheter and a soft closed-end trial catheter
- **available with Tip-in-Sight (TIS) technology for precise positioning of the catheter tip under ultrasound guidance**

LOADING CATHETER TIP OPTIONS

CLASSIC

TIP-IN-SIGHT

GYNO INFO SEMTRAC SET 5

Guiding catheter (A)	length	178 mm
	tube OD	2.55 mm

Loading catheter (C)	length	252 mm
	tube OD	1.70 mm
	total volume	0.22 ml
	total volume minus hub	0.089 ml

Sterile	ETO sterilized
---------	----------------


Packaging	75/shipping box
-----------	-----------------

Quality control	MEA and LAL
-----------------	-------------

Order numbers

- | | |
|------------------------------|------------------------|
| - Semtrac Set 5 | #2000 Set 5 |
| - Semtrac Set 5 Trial | #2000 Set 5 T |
| - Semtrac Set 5 Tip-in-Sight | #2000 Set 5 TIS |

BEFORE USING REFER TO INSTRUCTIONS FOR USE (IFU) FOR INDICATIONS, CONTRA-INDICATIONS AS WELL AS WARNINGS AND PRECAUTIONS.

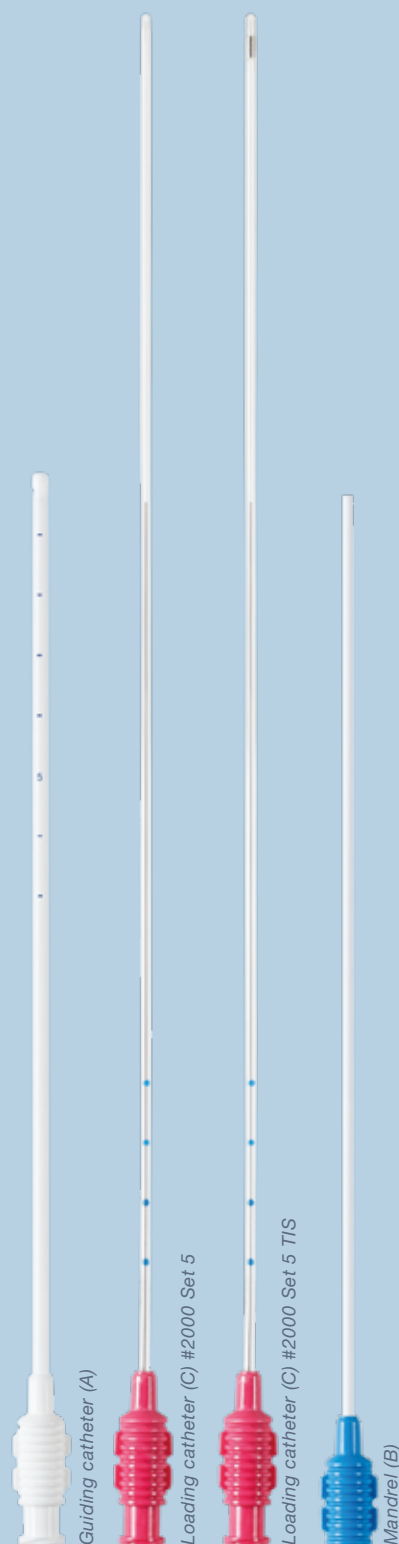

GMP Belgium Semtrac Set 5 06/2018



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be



Embryo Transfer

SEMTRAC SET W-180

SOFT EMBRYO TRANSFER CATHETER

The Semtrac Set is a soft embryo transfer catheter set for difficult embryo transfers.

FEATURES

- soft bullet shaped tip to minimize trauma
- 5 **free** mandrels (B), to preshape the guiding catheter, in every box of 25 pieces
- handle with luer lock for a more secure fitting to syringes
- trial catheter set available for cervix examination prior to embryo transfer consisting of a guiding catheter and a soft closed-end trial catheter
- **available with Tip-in-Sight (TIS) technology for precise positioning of the catheter tip under ultrasound guidance**


LOADING CATHETER TIP OPTIONS



GYNO INFO SEMTRAC SET W-180

Guiding catheter (A) length	141 mm
tube OD	2.55 mm
Loading catheter (C) length	214 mm
tube OD	1.70 mm
total volume	0.23 ml
total volume minus hub	0.074 ml
Sterile	ETO sterilized
Packaging	75/shipping box
Quality control	MEA and LAL
Order numbers	
- Semtrac Set W-180	#2000 W-180
- Semtrac Set W-180 Trial	#2000 W-180 T
- Semtrac Set W-180 Tip-in-Sight	#2000 W-180 TIS

BEFORE USING REFER TO INSTRUCTIONS FOR USE (IFU) FOR INDICATIONS, CONTRA-INDICATIONS AS WELL AS WARNINGS AND PRECAUTIONS.

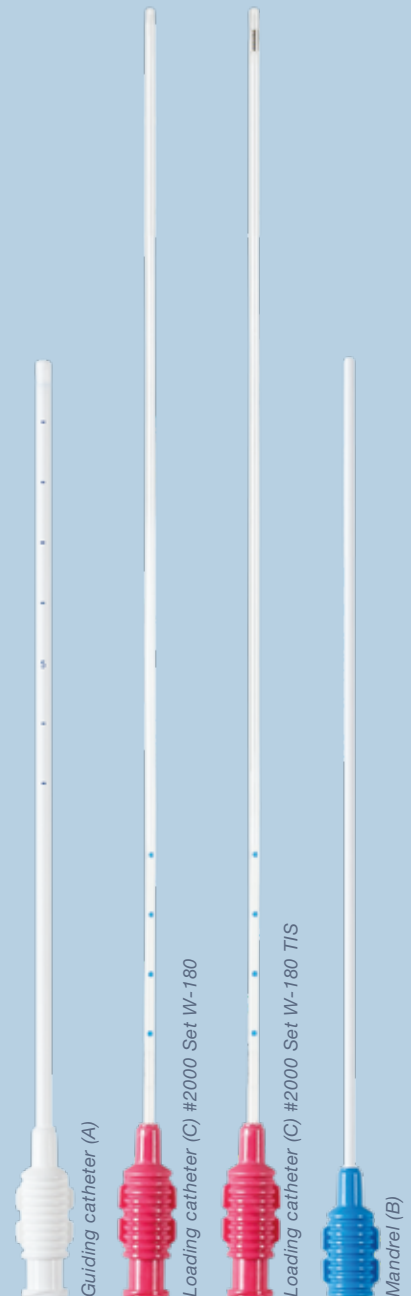

GMP Belgium Semtrac W-180 06/2018



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be



Embryo Transfer

SEMTRAC SET W-240

SOFT EMBRYO TRANSFER CATHETER

The Semtrac Set is a soft embryo transfer catheter set for difficult embryo transfers.

FEATURES

- soft bullet shaped tip to minimize trauma
- 5 **free** mandrels (B), to preshape the guiding catheter, in every box of 25 pieces
- handle with luer lock for a more secure fitting to syringes
- trial catheter set available for cervix examination prior to embryo transfer consisting of a guiding catheter and a soft closed-end trial catheter
- **available with Tip-in-Sight (TIS) technology for precise positioning of the catheter tip under ultrasound guidance**

LOADING CATHETER TIP OPTIONS


CLASSIC

TIP-IN-SIGHT

GYNO INFO SEMTRAC SET W-240

Guiding catheter (A) length	203 mm
tube OD	2.55 mm
Loading catheter (C) length	277 mm
tube OD	1.70 mm
total volume	0.26 ml
total volume minus hub	0.098 ml
Sterile	ETO sterilized
Packaging	75/shipping box
Quality control	MEA and LAL
Order numbers	
- Semtrac Set W-240	#2000 W-240
- Semtrac Set W-240 Trial	#2000 W-240 T
- Semtrac Set W-240 Tip-in-Sight	#2000 W-240 TIS

BEFORE USING REFER TO INSTRUCTIONS FOR USE (IFU) FOR INDICATIONS, CONTRA-INDICATIONS AS WELL AS WARNINGS AND PRECAUTIONS.

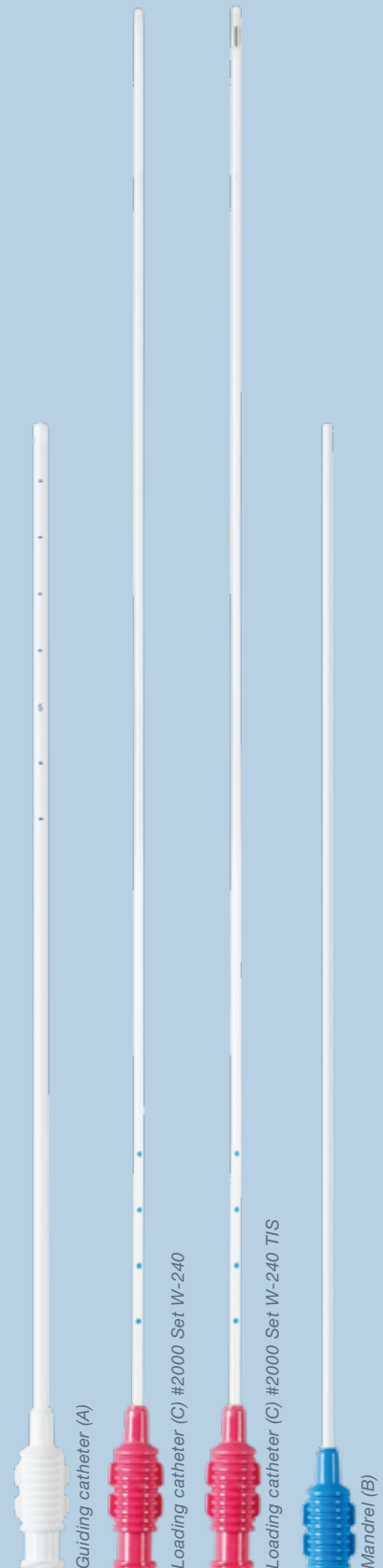

GMP Belgium Semtrac W-240 06/2018



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be



Embryo Transfer

THE TULIP

PRESHAPED EMBRYO TRANSFER CATHETER WITH TULIP SHAPED BULB

The Tulip is a set for difficult cervical passage, consisting of two parts: a guiding catheter (A) and a loading catheter (C).

FEATURES

- preshaped guiding catheter, ideally curved for easy entry when uterus is flexed
- the Tulip tip enables easy atraumatic entry through the cervical canal
- positioner for depth adjustment
- markings on outer tube for easy positioning
- handle with luer lock fitting for a more secure connection to syringes
- reinforced metal support of loading catheter for optimal control
- arrow marking on outer tube for curve indication

ORIGINAL PICTURE

GYNO INFO THE TULIP

Length
A) 200 mm
C) 266 mm

Tube OD
A) 2.20 mm
C) 1.00 mm

Total volume
Total volume minus hub
C) 0.21 ml
C) 0.069 ml


Tulip bulb
3 mm

Sterile
sterilized by radiation

Packaging
75/shipping box

Quality control
MEA and LAL

Order number
#4000 The Tulip

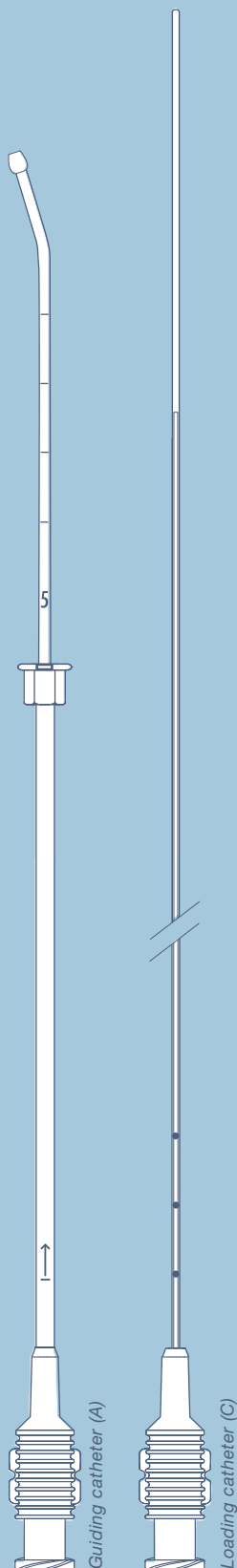

GMP Belgium The Tulip 04/2011



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS®
Practise better Medicine

THE TULIP MEMO

EMBRYO TRANSFER CATHETER WITH TULIP SHAPED BULB

The Tulip is a set for difficult cervical passage, consisting of two parts: a guiding catheter (A) and a loading catheter (C).

FEATURES

- the Tulip tip enables easy atraumatic entry through the cervical canal
- mandrel available for extra difficult transfers
- positioner for depth adjustment
- a combination of a soft, flexible intrauterine catheter and a solid cervix catheter with memory
- markings on outer tube for easy positioning
- memo function of the guiding catheter
- 5 free mandrels to preshape the guiding catheter in every box of 25 pieces
- handle with luer lock fitting for a more secure connection to syringes
- reinforced metal support of loading catheter for optimal control



ORIGINAL PICTURE

GYNO INFO THE TULIP MEMO

Length A) 200 mm
C) 266 mm

Tube OD A) 2.20 mm
C) 1.00 mm

Total volume C) 0.21 ml
Total volume minus hub C) 0.069 ml


Tulip bulb 3 mm

Sterile sterilized by radiation

Packaging 75/shipping box

Quality control MEA and LAL

Order number **#4000 Memo**

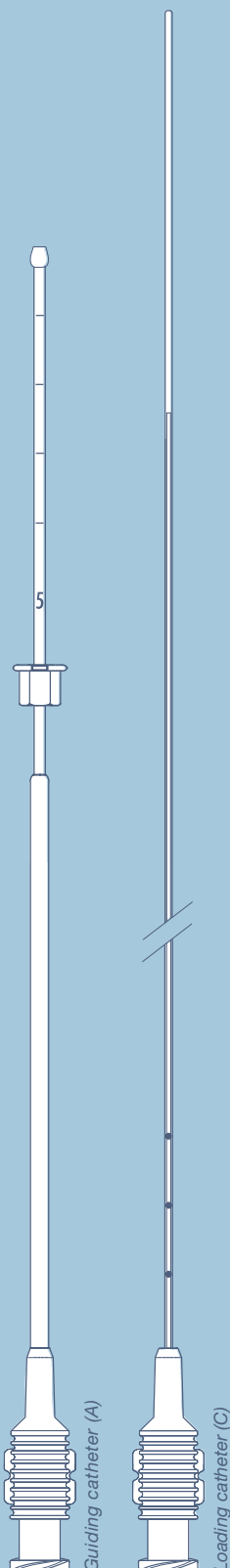

GMP Belgium The Tulip Memo 04/2011



EMBRYO TRANSFER CATHETER SETS WITH BULB SHAPED TIP

The renewed Tulip Embryo Transfer catheter range offers a variety of possibilities to support you enhancing your practice.

- The **bulb shaped** guiding catheters assist atraumatic entry into the cervical canal and are available in different shapes.
- The **Memo function** allows shaping the guiding catheter in any desired silhouette to suit various patients' anatomies.
- The Tulip **loading catheters** are available with a Classic tip or a Tip-in-Sight (TIS), which supports precise positioning under ultrasound guidance.


Guiding Catheter **bulb shape** options



The Tulip S
Oval shaped bulb



The Tulip L
Rigid round shaped bulb

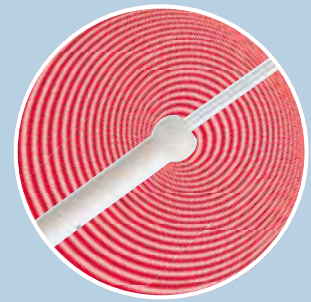

The Tulip LS
SOFT round shaped bulb

Loading Catheter **tip** options



Classic



Tip-in-Sight (TIS)
Echogenic




CHOOSE YOUR TULIP SET

PRODUCT		ORDER NUMBER	FEATURES			
			Echogenic tip	Preshaped with curve indicator	Memory function	Mandrel available 5 per 25 sets
 Oval shaped bulb	THE TULIP S	#4000 Tulip S		✓		
	THE TULIP S Memo	#4000 S Memo			✓	✓
	THE TULIP S Tip-in-Sight	#4000 Tulip S TIS	✓	✓		
	THE TULIP S Memo Tip-in-Sight	#4000 S Memo TIS	✓		✓	✓
 Rigid round shaped bulb	THE TULIP L	#4000 Tulip L		✓		
	THE TULIP L Memo	#4000 L Memo			✓	✓
	THE TULIP L Tip-in-Sight	#4000 Tulip L TIS	✓	✓		
	THE TULIP L Memo Tip-in-Sight	#4000 L Memo TIS	✓		✓	✓
 SOFT round shaped bulb	THE TULIP LS	#4000 Tulip LS		✓		
	THE TULIP LS Memo	#4000 LS Memo			✓	✓
	THE TULIP LS Tip-in-Sight	#4000 Tulip LS TIS	✓	✓		
	THE TULIP LS Memo Tip-in-Sight	#4000 LS Memo TIS	✓		✓	✓

GYNO INFO


Length	guiding catheter: 200 mm loading catheter: 266 mm
Tube OD	guiding catheter: 2.20 mm loading catheter: 1.00 mm
Total volume Total volume minus hub	loading catheter: 0.21 ml loading catheter: 0.069 ml
Bulb shape	oval shaped: 3 mm rigid/soft round shaped: 2.3 mm
Sterile	sterilized by gamma irradiation
Packaging	75/shipping box
Quality control	MEA and LAL

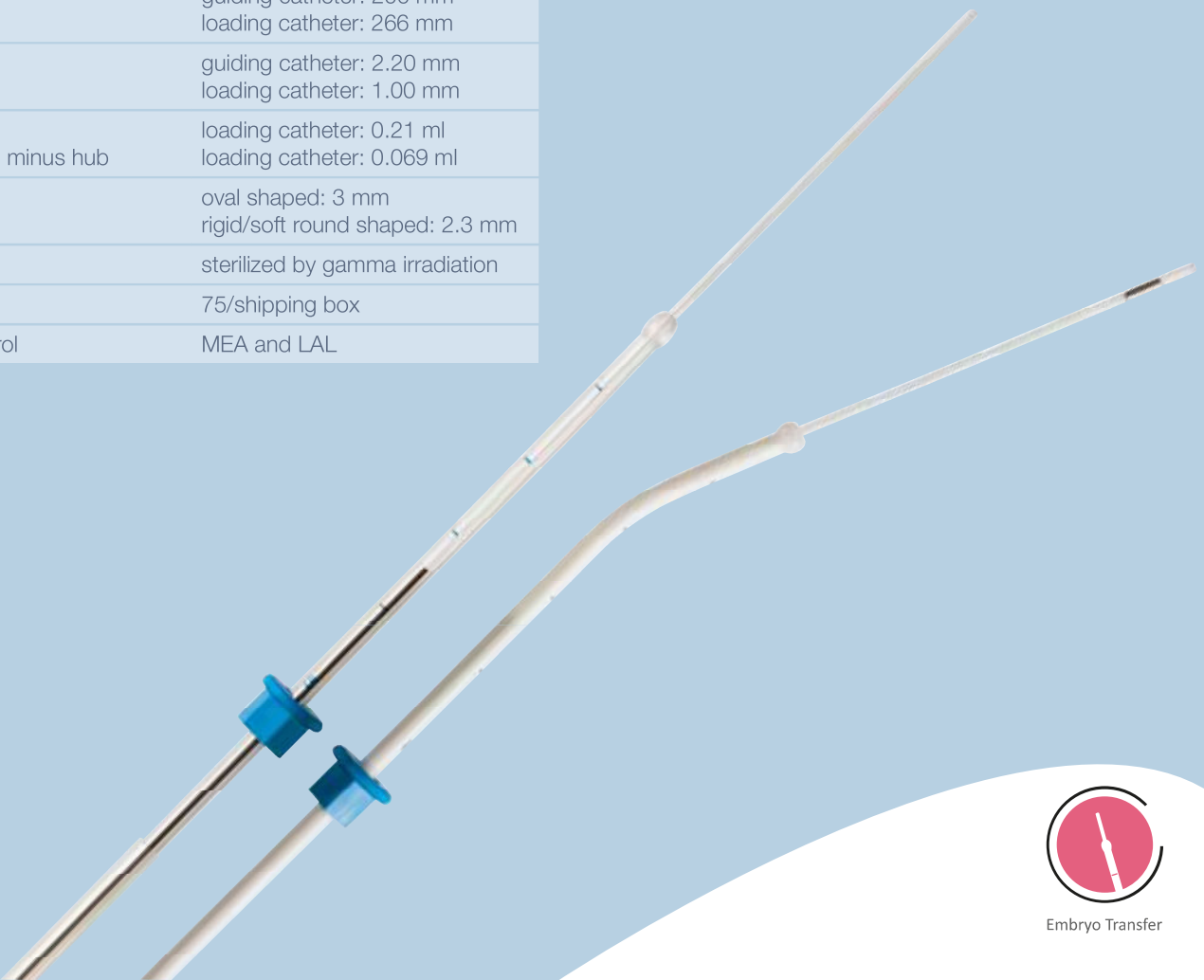

GYNETICS®

Practise better Medicine

FOLLICLE ASPIRATION NEEDLES

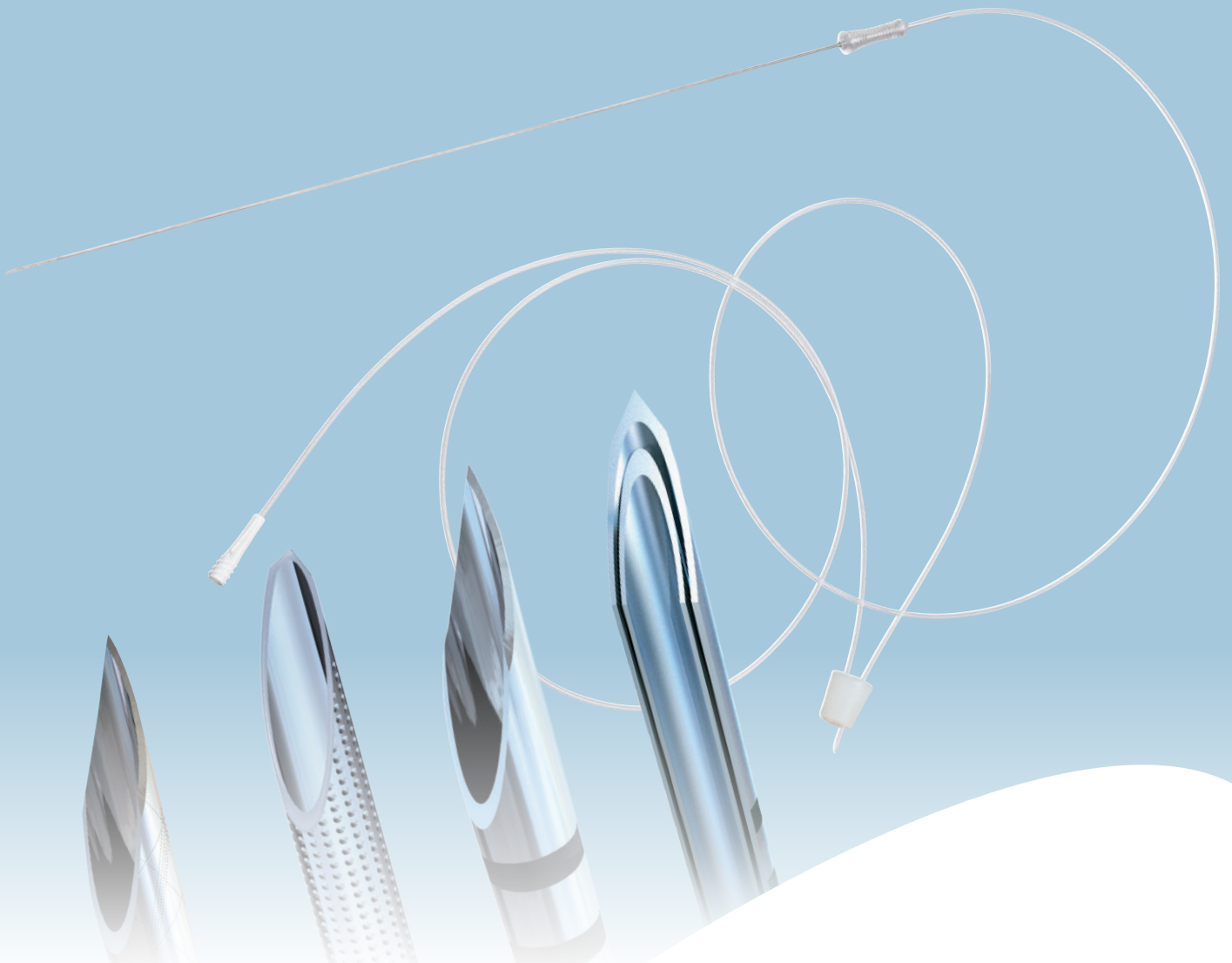



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine

GYNO INFO

SUPERIOR FOLLICLE ASPIRATION NEEDLE RANGE

PRODUCT	PAGE
Features and accessories	4
Single lumen	5-6
Single lumen flush	7-8
Double lumen flush	9-10



PREMIUM FOLLICLE ASPIRATION NEEDLE RANGE

PRODUCT	PAGE
Single lumen	11-12



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine

SUPERIOR FAS

FOLLICLE ASPIRATION SET

NEEDLE TYPES

3 needle types, available in various sizes and dimensions to determine your own specification:

Single Lumen

Single Lumen Flush

Double Lumen Flush

NEEDLE TIPS

2 types of needle tips to guarantee clean facet cutting:

Back Cut Bevel



Lancet Bevel



TIP MARKINGS

2 types of tip markings to offer excellent ultrasound visibility of the needle tip:

Dual Laser Marking




Embossed Marking




HANDLES

2 types of handles available for individual comfort (only for single lumen flush & double lumen flush):

U shaped handle

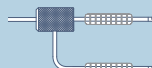

Y shaped handle

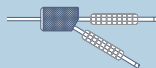

TRANSLUCENT TUBING

3 types of translucent tubing in any desired length (maximum length thick wall tubing is 120 cm / 47.2 inch):

Regular Tubing

for standard oocyte aspiration set up

Soft Tubing

allows control of the fluid

Thick Wall Tubing

no kinking of the tubing

CORKS

2 types of silicone corks for a secure connection to test tubes:

Small

14 x 18 x 20 mm

0.55 x 0.71 x 0.79 inch



Medium

21 x 27 x 30 mm

0.83 x 1.06 x 1.18 inch



Contact our sales department for other or special enquiries!



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine

SUPERIOR FAS SINGLE LUMEN

CUSTOM-MADE FOLLICLE ASPIRATION SET

The Superior FAS, the follicle aspiration set is implemented for ultrasound-guided, transvaginal harvesting of oocytes from follicles for IVF and ICSI.

FEATURES


- sharp lancet and back cut bevel tip available that ensure a clean facet cutting
- dual laser and embossed echogenic marked tip available which offer stable and excellent vision during ultrasound
- various needle sizes for ovaries with complicated and/or limited access
- translucent tubing in the same size as needle for safe oocyte retrieval
- arrow marking on handle for bevel point indication
- handle shape facilitates rotation
- custom-made



ORIGINAL PICTURE

GYNO INFO SUPERIOR FAS SINGLE LUMEN

Sterile	sterilized by gamma irradiation
Packaging	50/shipping box (standard configuration)
Quality control	MEA and LAL
Order number	#4551 (complete the configurator and receive your order number)

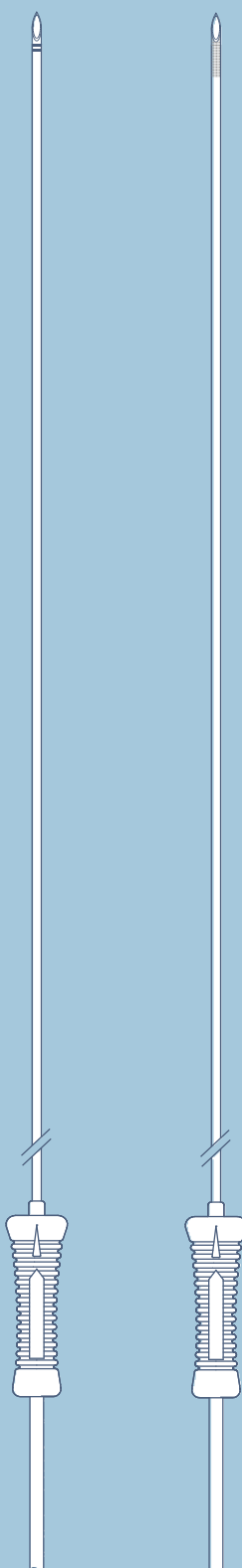


Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine

SUPERIOR FAS SINGLE LUMEN

CONFIGURE YOUR OWN CUSTOM-MADE SUPERIOR FOLLICLE ASPIRATION SET IN 6 STEPS

COPY OR DOWNLOAD THIS CONFIGURATOR AT WWW.GYNETICS.BE,
COMPLETE AND RETURN FOR INQUIRY TO RECEIVE YOUR CUSTOMIZED ORDER NUMBER

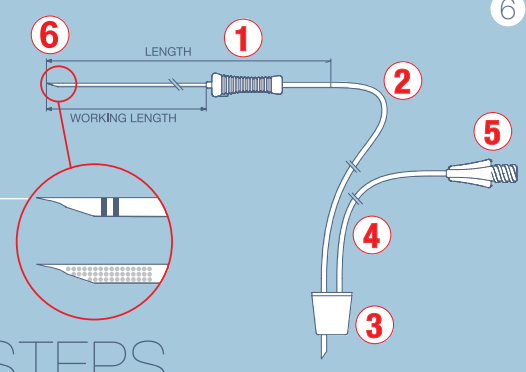

STEP 1 CHOOSE YOUR NEEDLE

TICK BOX ☒

Gauge	OD/ID (mm)	OD/ID (inch)	Length (mm)	Length (inch)	Control	
16G	1.6/1.3	0.063/0.051	400	15.75	Medium	☐
17G	1.5/0.9	0.059/0.035	400	15.75	Super rigid	☐
17G	1.5/0.9	0.059/0.035	350	13.78	Super rigid	☐
17G	1.4/1.0	0.055/0.039	400	15.75	Rigid	☐
17G	1.4/1.0	0.055/0.039	350	13.78	Rigid	☐
17G	1.4/1.2	0.055/0.047	400	15.75	Flexible	☐
17G	1.4/1.2	0.055/0.047	340	13.39	Flexible	☐
18G	1.2/1.0	0.047/0.039	300	15.75	Medium	☐
18G	1.2/1.0	0.047/0.039	350	13.78	Medium	☐

The working length of each needle is the length minus 50 mm/1.968 inch

STEP 2 CHOOSE YOUR TUBING FROM NEEDLE TO CORK

Type	Length needed	Remarks	
Regular Tubing	_____ cm/inch	n/a	☐
Soft Tubing	_____ cm/inch	n/a	☐
Thick Wall Tubing	_____ cm/inch	max. 120 cm/47.2 inch	☐

STEP 3 CHOOSE YOUR SILICONE CORK

Type	Size (mm)	Size (inch)	
Small Cork	14 x 18 x 20	0.55 x 0.71 x 0.79	☐
Medium Cork	21 x 27 x 30	0.83 x 1.06 x 1.18	☐

STEP 4 CHOOSE YOUR TUBING FROM CORK TO PUMP

Type	Length needed	Remarks	
Regular Tubing	_____ cm/inch	max 200 cm/ 78.7 inch	☐
Soft Tubing	_____ cm/inch	max 250 cm/ 98.4 inch	☐
Thick Wall Tubing	_____ cm/inch	max. 120 cm/47.2 inch	☐

YOUR CONNECTION TO THE PUMP

Female Luer Connector

STEP 5 CHOOSE YOUR VARIATION OF THE TIP

Back cut bevel with embossed echogenic marking		☐
Lancet bevel with embossed echogenic marking		☐
Back cut bevel with dual laser echogenic marking		☐
Lancet bevel with dual laser echogenic marking		☐

Contact our sales department for other or special enquiries!

SUPERIOR FAS SINGLE LUMEN FLUSH

CUSTOM-MADE FOLLICLE ASPIRATION SET

The Superior FAS, the follicle aspiration set is implemented for ultrasound-guided, transvaginal harvesting of oocytes from follicles for IVF and ICSI.


FEATURES

- flushing connector enables flushing and aspirating during complicated oocyte retrieval or when few mature oocytes occur
- flushing and aspirating procedure is performed separately through one lumen
- sharp lancet and back cut bevel tip available that ensure a clean facet cutting
- dual laser marked tip available which offer stable and excellent vision during ultrasound
- translucent tubing in the same size as needle for safe oocyte retrieval
- various needle sizes for ovaries with complicated and/or limited access
- handle shape facilitates rotation
- custom-made



GYNO INFO SUPERIOR FAS SINGLE LUMEN FLUSH

Sterile	sterilized by gamma irradiation
Packaging	50/shipping box (standard configuration)
Quality control	MEA and LAL
Order number	#4551 (complete the configurator and receive your order number)

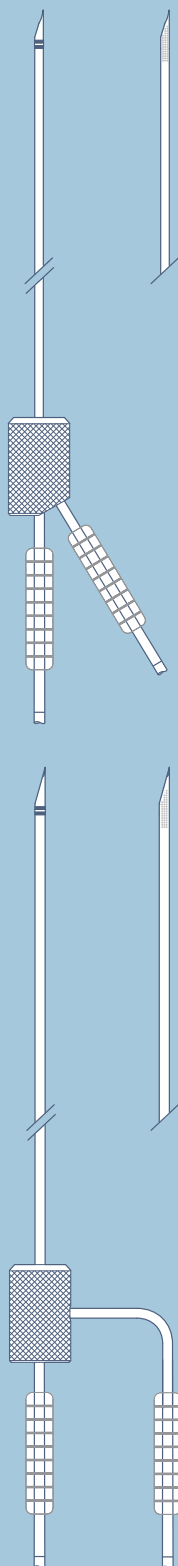


Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine

SUPERIOR FAS SINGLE LUMEN FLUSH



CONFIGURE YOUR OWN CUSTOM-MADE SUPERIOR FOLLICLE ASPIRATION SET IN 7 STEPS

COPY OR DOWNLOAD THIS CONFIGURATOR AT WWW.GYNETICS.BE,
COMPLETE AND RETURN FOR INQUIRY TO RECEIVE YOUR CUSTOMIZED ORDER NUMBER

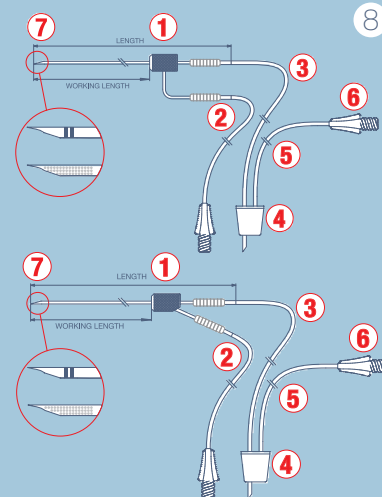

STEP 1 CHOOSE YOUR NEEDLE AND HANDLE

☒ TICK BOX

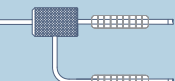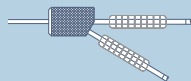Gauge	OD/ID (mm)	OD/ID (inch)	Length (mm)	Length (inch)	Control	
17G	1.40/1.20	0.055/0.047	352	13.858	Flexible	☐
17G	1.40/1.00	0.055/0.039	352	13.858	Rigid	☐
19G	2.20/1.62	0.087/0.063	352	13.858	Rigid	☐
The working length of each needle is 300 mm/11.811 inch						
U shaped handle		☐			Y shaped handle	 ☐

STEP 2 CHOOSE YOUR FLUSH TUBING

Tubing	Length needed	Remarks	
Regular Tubing	cm/inch	max. 195 cm/ 76.7 inch	☐
Soft Tubing	cm/inch	max. 160 cm/ 62.9 inch	☐
Thick Wall Tubing	cm/inch	max. 120 cm/47.2 inch	☐

STEP 3 CHOOSE YOUR TUBING FROM NEEDLE TO CORK

Tubing	Length needed	Remarks	
Regular Tubing	cm/inch	n/a	☐
Soft Tubing	cm/inch	n/a	☐
Thick Wall Tubing	cm/inch	max. 120 cm/47.2 inch	☐

STEP 4 CHOOSE YOUR SILICONE CORK

Type	Size (mm)	Size (inch)	
Small Cork	14 x 18 x 20	0.55 x 0.71 x 0.79	☐
Medium Cork	21 x 27 x 30	0.83 x 1.06 x 1.18	☐

STEP 5 CHOOSE YOUR TUBING FROM CORK TO PUMP

Type	Length needed	Remarks	
Regular Tubing	cm/inch	max. 200 cm/ 78.7 inch	☐
Soft Tubing	cm/inch	max. 250 cm/ 98.4 inch	☐
Thick Wall Tubing	cm/inch	max. 120 cm/47.2 inch	☐

YOUR CONNECTION TO THE PUMP

Female Luer Connector

STEP 6 CHOOSE YOUR VARIATION OF THE TIP

Back cut bevel with dual laser echogenic marking		☐
Lancet bevel with dual laser echogenic marking		☐

Contact our sales department for other or special enquiries!



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine

SUPERIOR FAS DOUBLE LUMEN FLUSH

CUSTOM-MADE FOLLICLE ASPIRATION SET

The Superior FAS, the follicle aspiration set is implemented for ultrasound-guided, transvaginal harvesting of oocytes from follicles for IVF and ICSI.

FEATURES



- flushing connector enables flushing and aspirating during complicated oocyte retrieval or when few mature oocytes occur
- double lumen allows simultaneous flushing and aspirating
- dual laser and embossed echogenic marked tip available which offer stable and excellent vision during ultrasound
- translucent tubing in the same size as needle for safe oocyte retrieval
- various needle sizes for ovaries with complicated and/or limited access
- handle shape facilitates rotation
- custom-made



ORIGINAL PICTURE

GYNO INFO SUPERIOR FAS DOUBLE LUMEN FLUSH

Sterile	sterilized by gamma irradiation
Packaging	50/shipping box (standard configuration)
Quality control	MEA and LAL
Order number	#4551 (complete the configurator and receive your order number)

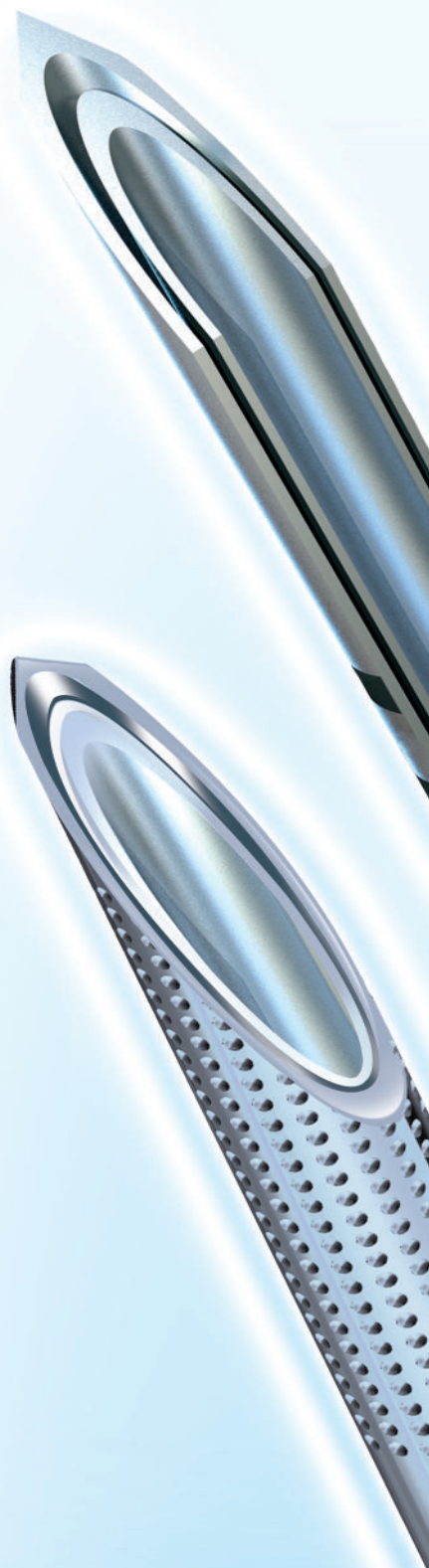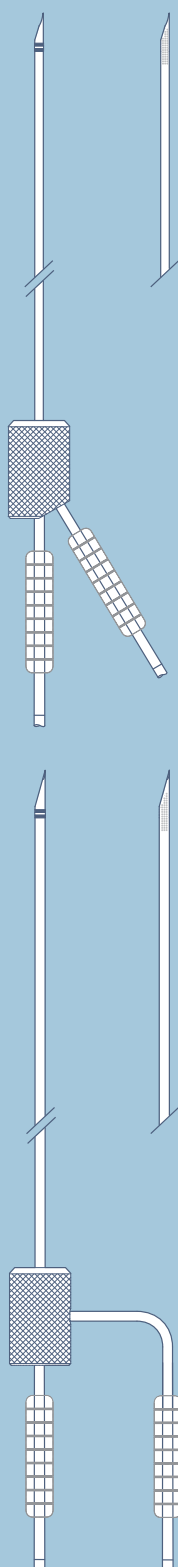


Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine

SUPERIOR FAS DOUBLE LUMEN FLUSH

CONFIGURE YOUR OWN CUSTOM-MADE SUPERIOR FOLLICLE ASPIRATION SET IN 7 STEPS

COPY OR DOWNLOAD THIS CONFIGURATOR AT WWW.GYNETICS.BE,
COMPLETE AND RETURN FOR INQUIRY TO RECEIVE YOUR CUSTOMIZED ORDER NUMBER


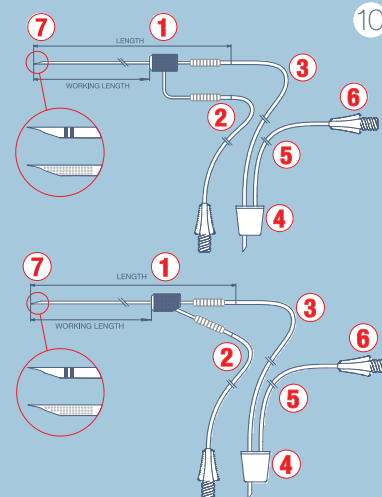


STEP 1 CHOOSE YOUR NEEDLE AND HANDLE

TICK BOX ☒

Gauge	OD/ID (mm)	OD/ID (inch)	Length (mm)	Length (inch)	Control	Handle	
16G	outer lumen 1.65/1.43 inner lumen 1.20/1.00	0.065/0.055 0.047/0.039	352	13.858	Rigid	U shaped	☐
17G	outer lumen 1.50/1.20 inner lumen 1.00/0.80	0.059/0.047 0.039/0.032	352	13.858	Rigid	Y shaped	☐

The working length of each needle is 300 mm/11.811 inch

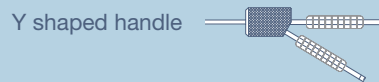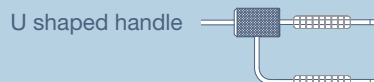

STEP 2 CHOOSE YOUR FLUSH TUBING

Tubing	Length needed	Remarks	
Regular Tubing	_____ cm/inch	max. 195 cm/ 76.7 inch	☐
Soft Tubing	_____ cm/inch	max. 160 cm/ 62.9 inch	☐
Thick Wall Tubing	_____ cm/inch	max. 120 cm/47.2 inch	☐

STEP 3 CHOOSE YOUR TUBING FROM NEEDLE TO CORK

Tubing	Length needed	Remarks	
Regular Tubing	_____ cm/inch	n/a	☐
Soft Tubing	_____ cm/inch	n/a	☐
Thick Wall Tubing	_____ cm/inch	max. 120 cm/47.2 inch	☐

STEP 4 CHOOSE YOUR SILICONE CORK

Type	Size (mm)	Size (inch)	
Small Cork	14 x 18 x 20	0.55 x 0.71 x 0.79	☐
Medium Cork	21 x 27 x 30	0.83 x 1.06 x 1.18	☐

STEP 5 CHOOSE YOUR TUBING FROM CORK TO PUMP

Type	Length needed	Remarks	
Regular Tubing	_____ cm/inch	max. 200 cm/ 78.7 inch	☐
Soft Tubing	_____ cm/inch	max. 250 cm/ 98.4 inch	☐
Thick Wall Tubing	_____ cm/inch	max. 120 cm/47.2 inch	☐

YOUR CONNECTION TO THE PUMP

Female Luer Connector

STEP 6 CHOOSE YOUR VARIATION OF THE TIP

Dual laser echogenic marking		☐
Embossed echogenic marking		☐

Contact our sales department for other or special enquiries!



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine

PREMIUM FAS SINGLE LUMEN

FOLLICLE ASPIRATION SET

The Premium FAS, the follicle aspiration set is implemented for ultrasound-guided, transvaginal harvesting of oocytes from follicles for IVF and ICSI.

FEATURES


- standard follicle aspiration needle with regular tubing
- sharp lancet and back cut bevel tip available that ensure a clean facet cutting
- diamond laser ultrasound markings on tip which offer stable and excellent vision during ultrasound
- translucent tubing in the same size as needle for safe oocyte retrieval
- arrow marking on handle for bevel point indication
- handle shape facilitates rotation



ORIGINAL PICTURE

GYNO INFO PREMIUM FAS SINGLE LUMEN

Sterile	sterilized by gamma irradiation
Packaging	50/shipping box (standard configuration)
Quality control	MEA and LAL
Order number	#4551 NS-AS1 #4551 NS-BS1 #4551 NS-CS1 #4551 NS-DS1 #4551 NS-AS1BCB #4551 NS-CS1BCB #4551 NS-DS1BCB

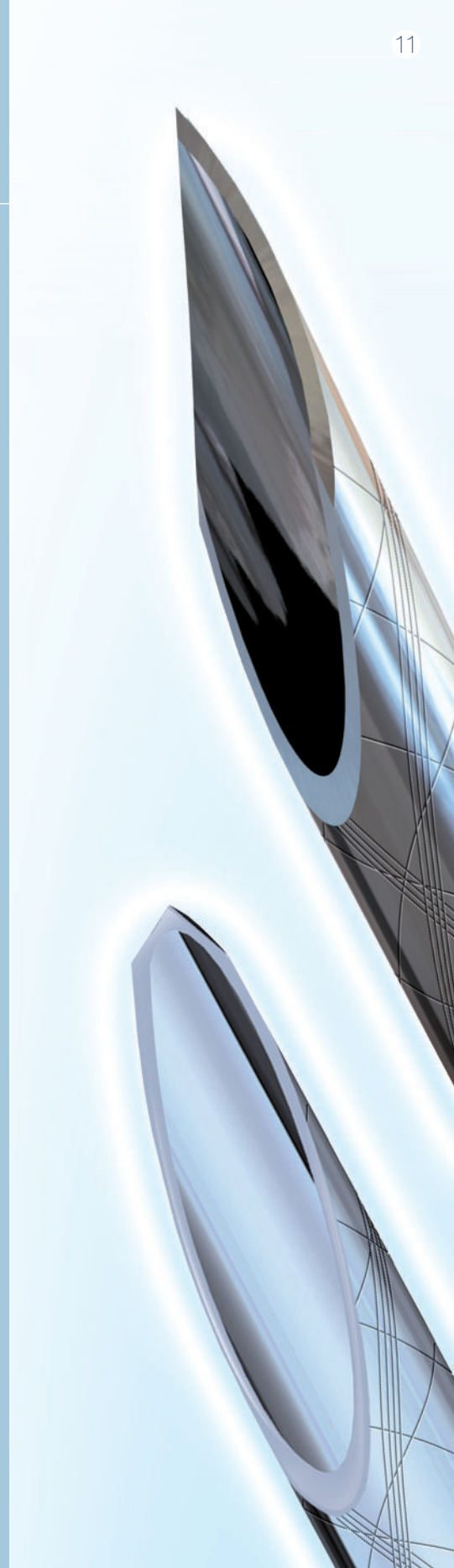

Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine


GYNO INFO

PREMIUM FAS SINGLE LUMEN

① Order number	Gauge	OD/ID (mm)	Length (mm)	Working length (mm)	Bevel	Control
#4551 NS-AS1	17G	1.4/1.2	400	350	Lancet	Flexible
#4551 NS-BS1	17G	1.4/1.2	340	290	Lancet	Flexible
#4551 NS-CS1	17G	1.4/1.0	400	350	Lancet	Rigid
#4551 NS-DS1	17G	1.4/1.0	340	290	Lancet	Rigid
#4551 NS-AS1BCB	17G	1.4/1.2	400	350	Back cut bevel	Flexible
#4551 NS-CS1BCB	17G	1.4/1.0	400	350	Back cut bevel	Rigid
#4551 NS-DS1BCB	17G	1.4/1.0	340	290	Back cut bevel	Rigid

Order number	Gauge	OD/ID (inch)	Length (inch)	Working length (inch)	Bevel	Control
#4551 NS-AS1	17G	0.055/0.047	15.75	13.779	Lancet	Flexible
#4551 NS-BS1	17G	0.055/0.047	13.39	11.417	Lancet	Flexible
#4551 NS-CS1	17G	0.055/0.039	15.75	13.779	Lancet	Rigid
#4551 NS-DS1	17G	0.055/0.039	13.39	11.417	Lancet	Rigid
#4551 NS-AS1BCB	17G	0.055/0.047	15.75	13.779	Back cut bevel	Flexible
#4551 NS-CS1BCB	17G	0.055/0.039	15.75	13.779	Back cut bevel	Rigid
#4551 NS-DS1BCB	17G	0.055/0.039	13.39	11.417	Back cut bevel	Rigid

- ② Tubing from needle to cork : Regular tubing 100 cm / 39.369 inch
- ③ Size silicone cork : Small 14 x 18 x 20 mm / 0.551 x 0.709 x 0.787 inch
- ④ Tubing from cork to pump : Regular tubing 100 cm / 39.369 inch

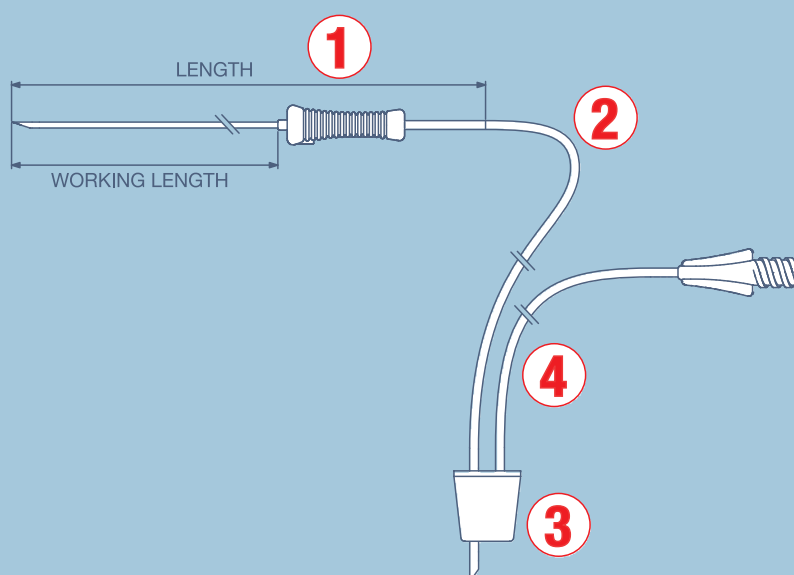

GYNO INFO

SAMPLING RANGE

PRODUCT	ORDER NUMBER
Detecta	#4163
Probet	#4164
Gynosampler	#9164
In/out	#4165

GMP Belgium Sampling 11/2017



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine


DETECTA

CERVICAL SCREENING DEVICE DESIGNED TO OBTAIN OPTIMAL SAMPLES

The Detecta is a simple but efficient cervical sampling device.

FEATURES

- allows in-office sampling
- systematic brushing of the transitional zone
- helps to detect neoplasia at an early stage
- easy to handle

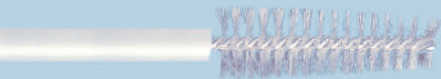

ORIGINAL PICTURE

GYNO INFO DETECTA

Length	194 mm
Packaging	1250/shipping box
Order number	#4163



GMP Belgium Detecta 06/09



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine

PROBET

ENDOMETRIAL HISTOLOGY

The Probet is an endometrial suction curette for a histologic lining biopsy.

FEATURES

- double headed piston secures a steady vacuum
- tissue collection through side hole
- smooth rounded distal tip and small size for minimal patient discomfort
- depth markings for ease in locating depth of insertion in the uterus



ORIGINAL PICTURE

GYNO INFO PROBET

Length	262 mm
Tube OD	2.93 mm
Sterile	ETO sterilized
Packaging	75/shipping box
Order number	#4164



GMP Belgium Probet 04/2011



GYNOSAMPLER

ENDOMETRIAL HISTOLOGY

The Gynosampler is the American name for the Probet which is sold in the USA.

The Gynosampler is an endometrial suction curette for a histologic lining biopsy.

FEATURES

- double headed piston secures a steady vacuum
- tissue collection through side hole
- smooth rounded distal tip and small size for minimal patient discomfort
- depth markings for ease in locating depth of insertion in the uterus



ORIGINAL PICTURE

GYNO INFO GYNOSAMPLER

Length	262 mm
Tube OD	2.93 mm
Sterile	ETO sterilized
Packaging	250/shipping box (US only)
Order number	#9164 (US only)



GMP Belgium Gynosampler 04/2017



IN / OUT

ENDOMETRIAL HISTOLOGY

The IN/OUT is an endometrial suction curette for a histologic lining biopsy.

FEATURES

- double headed piston secures a steady vacuum
- tissue collection through the tip
- smooth rounded distal tip and small size for minimal patient discomfort
- depth markings for ease in locating depth of insertion in the uterus



ORIGINAL PICTURE

GYNO INFO IN/OUT

Length	262 mm
Tube OD	2.93 mm
Sterile	ETO sterilized
Packaging	75/shipping box
Order number	#4165



GMP Belgium IN/OUT 04/2011



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine

GYNO INFO

IVF LABORATORY RANGE

PRODUCT

ORDER NUMBER

Handling Pipette

#3145
#3100
#3125
#3145 TP
#3100 TP
#3125 TP

Strip-ID

#3000 Strip-ID125/10
#3000 Strip-ID130/10
#3000 Strip-ID135/10
#3000 Strip-ID140/10
#3000 Strip-ID145/10
#3000 Strip-ID150/10
#3000 Strip-ID170/10
#3000 Strip-ID175/10
#3000 Strip-ID200/10
#3000 Strip-ID275/10
#3000 Strip-ID300/10

Pipet-ID

#3000 Pipet-ID
#3000 Pipet-IDs

GMP Belgium IVF Laboratory 09/2016



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine

HANDLING PIPETTE

PIPETTE FOR HANDLING EMBRYOS AND OOCYTES

The handling pipette is a single-use IVF laboratory micropipette used for transferring oocytes and embryos from one dish to another.

FEATURES


- IVF laboratory micropipette for dish-to-dish transfer of oocytes and embryos
- will not scratch dishes
- available in three lengths: 45, 55 and 94 mm
- available in two packaging variations:
single packed and packed by ten pieces



ORIGINAL PICTURE

GYNO INFO HANDLING PIPETTE

Length	45, 55 or 94 mm
Tube OD	1.25 mm
Tube ID	1.00 mm
Sterile	sterilized by gamma irradiation
Packaging	250/shipping box (1/pouch) 1250/shipping box (10/pouch)
Quality control	MEA and LAL
Order number	#3145 (45 mm, 1/pouch)) #3100 (55 mm, 1/pouch) #3125 (94 mm, 1/pouch) #3145 TP (45 mm, 10/pouch) #3100 TP (55 mm, 10/pouch) #3125 TP (94 mm, 10/pouch)

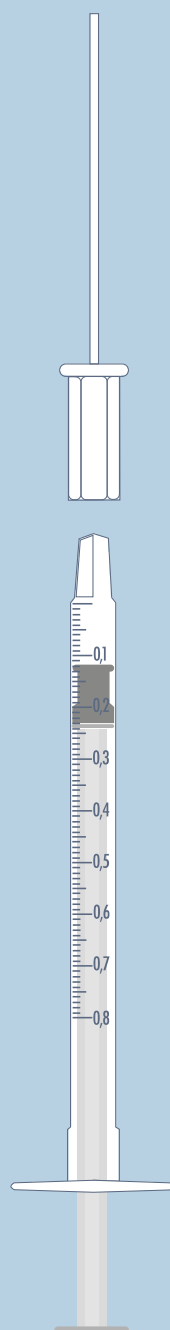

GMP Belgium Handling Pipette 09/2016



STRIP-ID

PIPETTE FOR OOCYTES MANIPULATION

The Strip-ID is a single-use IVF laboratory micropipette used for oocytes manipulation and handling of embryos and oocytes

FEATURES

- superior quality that ensures smooth edges for precise denudation
- strict individual microscopic quality inspection
- flexible, unbreakable material that will not scratch dishes
- convenient colour-coded multi-packaging that allows easy release of a single pipette
- extensive variety of sizes available to suit your application
- fits the "Pipet-ID(s)" and is compatible with most popular used micropipettors

ORIGINAL PICTURE



ORIG. PIC. 200X
ENLARGED

GYNO INFO STRIP-ID

Length	90 mm		
Inner diameter μm	125 μm	145 μm	200 μm
	130 μm	150 μm	275 μm
	135 μm	170 μm	300 μm
	140 μm	175 μm	
Sterile	sterilized by gamma irradiation		
Packaging	10 pipettes / vial 15 vials / inner box 75 vials of 10 pipettes / shipping box		
Quality control	MEA and LAL		
Order number	#3000 Strip-ID125/10 #3000 Strip-ID130/10 #3000 Strip-ID135/10 #3000 Strip-ID140/10 #3000 Strip-ID145/10 #3000 Strip-ID150/10 #3000 Strip-ID170/10 #3000 Strip-ID175/10 #3000 Strip-ID200/10 #3000 Strip-ID275/10 #3000 Strip-ID300/10		

GMP Belgium Strip-ID 12/2016



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine

GYNO INFO

STRIP-ID

Order number	ID	Colour
#3000 Strip-ID125/10	125 µm	blue
#3000 Strip-ID130/10	130 µm	yellow
#3000 Strip-ID135/10	135 µm	yellow
#3000 Strip-ID140/10	140 µm	white
#3000 Strip-ID145/10	145 µm	white
#3000 Strip-ID150/10	150 µm	black
#3000 Strip-ID170/10	170 µm	red
#3000 Strip-ID175/10	175 µm	red
#3000 Strip-ID200/10	200 µm	dark red
#3000 Strip-ID275/10	275 µm	green
#3000 Strip-ID300/10	300 µm	transparant



Easy release of a single pipette



Colour-coded multi-packaging



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine

PIPET-ID

MICROPIPETTOR

The Pipet-ID and Pipet-IDs are reusable IVF laboratory micropipettors developed to fit the "Strip-ID" micropipettes and most other popular micropipettes (Ø 0.9 mm OD)

FEATURES

- Available in two different sizes to fit every size of hand
 - **Pipet-ID** for regular to large hands
 - **Pipet-IDs** for small to regular hands
- Ergonomic design
- Adjustable aspiration volume (min ≈ 0.2 µl, max ≈ 3 µl)
- Total expulsion feature
- Stainless steel table holder included
- Spare parts kit included (includes extra seals, plunger wire and hexagon key)
- Designed for use with the "Strip-ID" micropipettes. Also fits most plastic micropipettes (Ø 0.9 mm OD, 0.6 mm ID)




ORIGINAL PICTURE

GYNO INFO PIPET-ID

Length (setting on position 3)	143 mm
Weight	57 gr
Order number	#3000 Pipet-ID

GYNO INFO PIPET-IDs

Length (setting on position 3)	120 mm
Weight	43 gr
Order number	#3000 Pipet-IDs

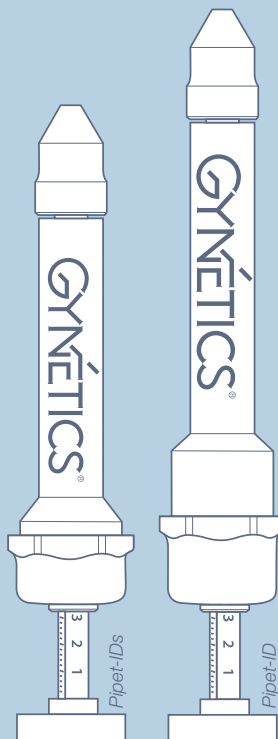

GMP Belgium Pipet-ID 10/2017



Colour differences may occur between batches due to variations in the colouring process and be considered as normal. GMP will not accept any claim or remark concerning these eventual differences.



Gynétics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNÉTICS®
Practise better Medicine