

THE TULIP

PRESHAPED EMBRYO TRANSFER CATHETER WITH TULIP SHAPED BULB

The Tulip is a set for difficult cervical passage, consisting of two parts: a guiding catheter (A) and a loading catheter (C).

FEATURES

- preshaped guiding catheter, ideally curved for easy entry when uterus is flexed
- the Tulip tip enables easy atraumatic entry through the cervical canal
- positioner for depth adjustment
- markings on outer tube for easy positioning
- handle with luer lock fitting for a more secure connection to syringes
- reinforced metal support of loading catheter for optimal control
- arrow marking on outer tube for curve indication

ORIGINAL PICTURE

GYNO INFO THE TULIP

Length A) 200 mm
C) 266 mm

Tube OD A) 2.20 mm
C) 1.00 mm

Total volume C) 0.21 ml
Total volume minus hub C) 0.069 ml


Tulip bulb 3 mm

Sterile sterilized by radiation

Packaging 75/shipping box

Quality control MEA and LAL

Order number #4000 The Tulip

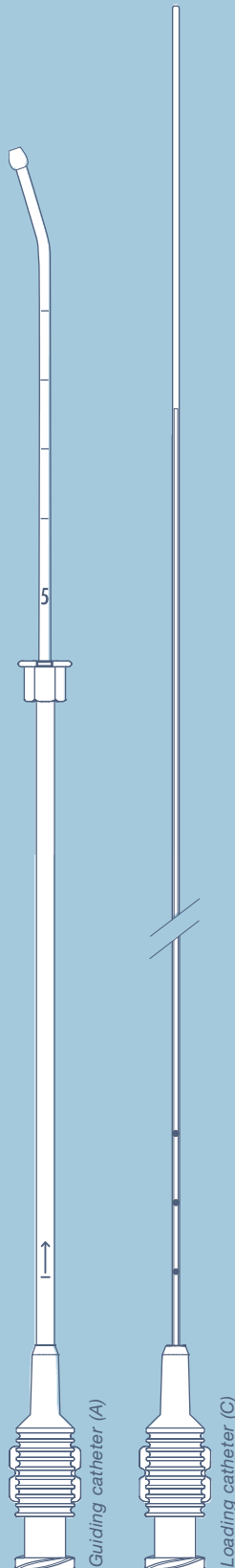

GMP Belgium The Tulip 04/2011



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine