#3000 PIPET-ID/ #3000 PIPET-IDS

Important information, read carefully before use

Manufactured by:
Gynétics Medical Products N.V.
Rembert Dodoensstraat 51
3920 Lommel, Belgium

HOW SUPPLIED

REF #3000 Pipet ID(s)
REF #3000 Spare Parts Kit(s)

DESCRIPTION

The Pipet-ID(s) is a non-sterile reusable micropipettor to fit the Strip ID micropipettes and most popular micropipettes (Ø .9 mm OD). Available in two different sizes to fit different size of hand.

INTENDED USE:

The Pipet-ID(s) is a reusable IVF laboratory micropipettor.

HOW TO USE

PREPARATION FOR USE

- Before first use, it is recommended to clean the Pipet-ID(s) by wiping with isopropyl alcohol or similar cleaning agents.
- The Pipet-ID(s) is suitable for ultrasonic[1,3] and autoclave cleaning[2,3].
 Please note: remove the 3 seals (injection cleaner, injection holder for the needle and injection clamp for the needle) before autoclave process.
- If maintenance is carried out properly, there is no restriction on the number of times the Pipet-ID(s) may be reused.

[1]50°C / 20 Minutes / Rinse with distilled water.
[2]Moist Heat / 125°C / 10 Minutes / 0.1 Mpa
[3]After ultrasonic or autoclave cycle let dry before use.

APPLYING A PIPETTE

- Before proceeding, make sure the collet is loosened one turn from the stop **(Off) (See A)**.
- Depress the plunger button until the injection needle appears and is long enough for applying a Pipette. **(See C)**
- Apply the Pipette over the injection needle and into the collet until it stops against the seals. Push the Pipette through the first resistance until it can no longer be pushed any further **(See A)** i.e. up to the internal stop.
- Turn the collet gently clock wards towards its stop **(On)** in order to fasten the pipette prior to denuding **(See A)**.
 Please note: Do not over-tighten the collet. This may affect the smoothness of the plunger during denudation.
- Check the fixation of the pipette before use by depressing the plunger button several times.

NORMAL USE

- Depress the plunger button to the first stop and hold it in this position (black ring stays visible) **(See B & C)**.
- Immerse the tip of the Pipette into the sample media while keeping the plunger depressed **(See C)**.
- Slowly release the plunger to aspirate the required amount of medium and sample **(See C)**.
- Expel the sample by completely depressing the plunger until there is no more media left in the Pipette **(See C)** or use the total expulsion feature **(See D)**.

REMOVING A PIPETTE

- Before removing a Pipette, check whether the sample has been expelled safely.
- To expel any residual media in the Pipette, the total expulsion feature can be used **(See D)**.
- Ensure the collet is loosened one turn from the stop **(Off)** before continuing **(See A)**.
- Depress the plunger until there appears at least 10 mm of injection needle inside the Pipette **(See B)**.
- Hold the Pipette firmly between your fingers and pull it gently out of the collet, using the injection needle as a guide **(See A)**.

TOTAL EXPULSION FEATURE

- The Pipet-ID(s) has a total expulsion capacity of 1.4 µl on top of the aspiration setting.
- To remove the last remaining sample media from the Pipette, depress the plunger until it cannot be depressed any further, and release **(See D)**.
- The total expulsion feature will be initiated if the black ring is pressed by the plunger. During normal use, the black ring should be avoided unless one deliberately wants to use the total expulsion feature.

ADJUSTING VOLUME

- The Pipet-ID(s) is set so that it can draw up in total 3 µl* of media.
- When you rotate the thumb rest clockwise the setting decreases (anti-clockwise increases).
- The settings are in steps of 0.2 µl*, marked with 3*, 2* and 1 µl* **(See B & E)**.
- To reset the Pipet-ID(s) to the maximum setting, the thumb rest should be turned anti-clockwise until the internal stop is reached **(See B & E)**.

**Indication only*

MAINTENANCE

REPLACING THE SEALS

The seals need to be replaced every six months or earlier if:

- The Pipette fills with media before the plunger is released.
- The Pipette still fills with media even though the collet was fastened.
- The seals look worn.

A replacement kit, including seals and a injection needle, is supplied with the Pipet-ID(s) (#3000 Spare Parts Kit(s)).

HOW TO REPLACE THE SEALS

- Take the Pipette out of the Pipet-ID(s) **(See A)**.
- Loosen the collet **(See F)**.
- Remove the old seals from the injection needle, or from the inside of the collet if they are stuck (the blunt end of a used Pipette can be used for this), and discard them.
- Inspect the injection needle and the inside of the collet, and clean and disinfect, if necessary.
- Slide the seals over the injection needle in the following order **(See F)**:
 1st Largest red seal (injection cleaner).
 2nd Translucent seal, largest hole facing away from the body of the Pipet-ID(s) (injection holder for the needle).
 3rd Smallest red seal (injection clamp for the needle).

REF 00499h
Last updated: 2026-07-02

- If the seals are applied in any other quantity, order or orientation as described above the Pipet-ID(s) will not function properly.

REPLACING THE INJECTION NEEDLE

The injection needle needs to be replaced if:

- It is not possible anymore to apply a Pipette over the injection needle.
- The injection needle is bent or damaged.
- It is not possible anymore to depress the plunger.

INSERTING A NEW INJECTION NEEDLE

- Take the Pipette out of the Pipet-ID(s) **(See A)**.
- Loosen the collet **(See F)**.
- Take out the three seals **(See F)**.
- Loosen the broader end of the Pipet-ID(s) body anti-clockwise until it can move over the plunger **(See G 1st step)**.
- Turn the thumb rest anti-clockwise to loosen the mechanism from the body **(See G 2nd step)**.
- Carefully loosen the two grub screws by using the hexagon key, supplied with the Pipet-ID(s), until the injection needle is released **(See G 3rd step)**.
- Take out the injection needle **(See G 4th step)**.
- Insert the new injection needle in the narrow opening in the mechanism and carefully push it in. Make sure that at least 60 mm remains visible.
- Fasten the two grub screws to tightly attach the new injection needle **(See G 3rd step)**.
- Gently insert the sliding nut back into the body of the Pipet-ID(s) by turning the thumb rest clockwise **(See G 2nd & 1st steps)**.
- Turn the thumb rest clockwise until the volume markings on the injection needle are no longer visible
- Bring the broader end of the Pipet-ID(s) body clockwise back into place **(See G)**.
- Turn the thumb rest anti-clockwise until it stops to reset the aspirate volume setting to its maximum of 3µl **(See E)**.
- Place back the three seals and the collet **(See F)** (1st injection cleaner, 2nd injection holder for the needle and 3rd injection clamp for the needle).
- Test the Pipet-ID(s) to guarantee a smooth movement and repeat steps if necessary.

CAUTION

- The denuding pipettes are for single use only
- The Pipet-ID(s) should only be used by qualified personnel.
- Before use, make sure the Pipet-ID(s) is properly cleaned and disinfected by wiping with isopropyl alcohol or similar cleaning agents.

TROUBLESHOOTING		
PROBLEM	**POTENTIAL CAUSE**	**SOLUTION**
Pipette does not fit	Incorrect Pipette size used	Only use plastic Pipettes with an outer diameter of 0.9 mm
	Seals not correctly fitted	Follow instructions 'Replacing the seals'
	Plunger not depressed	Follow instructions 'Normal use'
	Injection needle damaged	Replace injection needle
	Pipette damaged	Take another Pipette
Pipette falls out	Pipette is not fitted correctly	Follow instructions 'Applying a pipette' With attention to underlined part.
Reduced control of sample	Pipette not correctly fitted	Follow instructions 'Applying a Pipette'
	Pipette damaged	Take another Pipette
	Seals damaged	Replace seals
	Seals are not fitted in the right order or orientation	Follow instructions 'replace seals' with attention to underlined part
	Injection needle damaged	Replace injection needle
	Injection needle loose	Secure injection needle
	Sample delivery set to zero	Reset sample delivery to 1- 3 µl
	Collet too loose	Fasten collet
Plunger does not depress	Thumb rest tightened to body	Reset sample delivery to 1- 3 µl
	Injection needle damaged	Replace injection needle
	Seals damaged	Replace seals
	Injection needle loose	Secure injection needle
	Pipette not correctly fitted	Follow instructions 'Applying a Pipette'
Aspiration volume cannot be adjusted	Maximum or minimum reached	Follow instructions 'Adjusting volume'
	Not correctly maintained	Follow instructions 'Replacing the injection needle'
	Pipet-ID(s) damaged	Contact Gynétics

REF Catalogue number

LOT Batch code

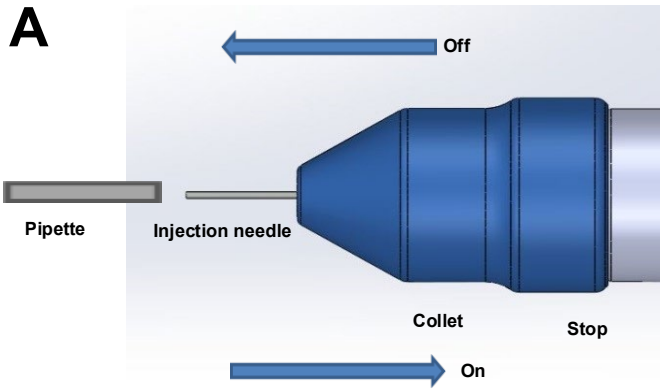
Manufacturer

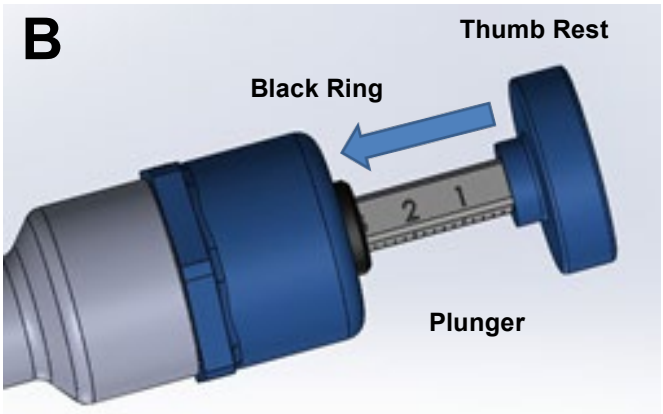
Consult instructions for use

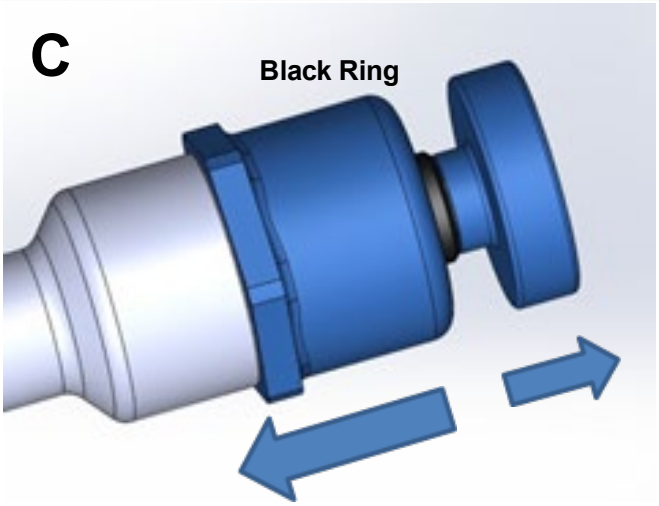
CE marking of conformity

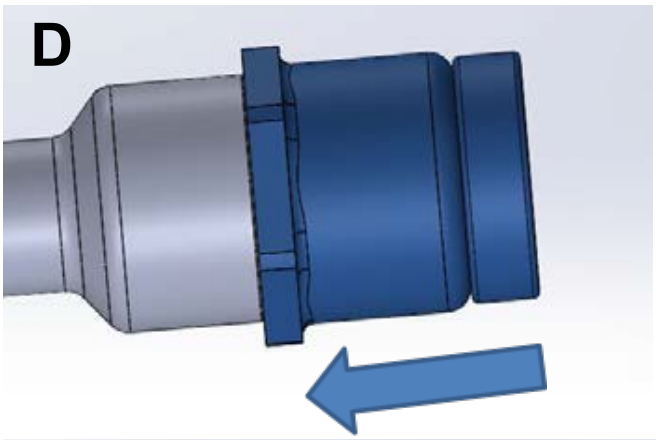
MD Medical device

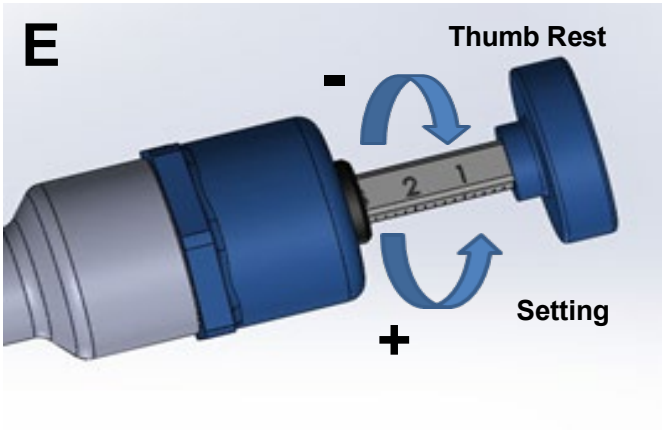
Keep dry

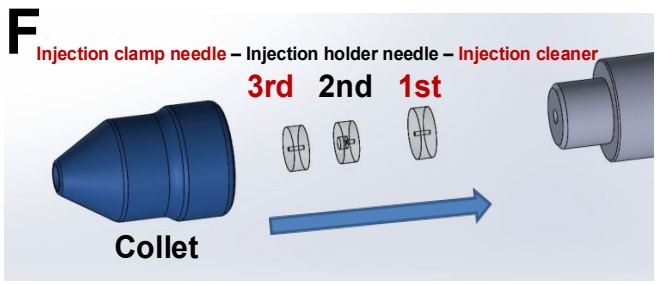
A
Off
Pipette
Injection needle
Collet
Stop
On

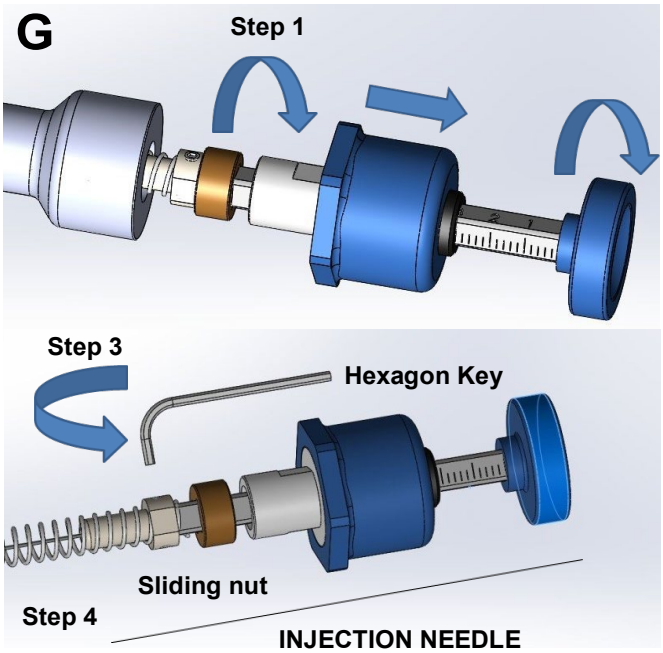
B
Thumb Rest
Black Ring
2
1
Plunger

C
Black Ring

D

E
Thumb Rest
-
2
1
Setting
+

F
Injection clamp needle – Injection holder needle – Injection cleaner
3rd
2nd
1st
Collet

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